CEPILLO DENTAL ELÉCTRICO SPAZZOLINO ELETRICO NBZ 45 A1

CEPILLO DENTAL ELÉCTRICO

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SPAZZOLINO ELETRICO

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ELECTRIC TOOTHBRUSH

Operating instructions

ELEKTRISCHE ZAHNBÜRSTE

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IAN 373276_2104

IAN 373276_2104

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
- 14 -

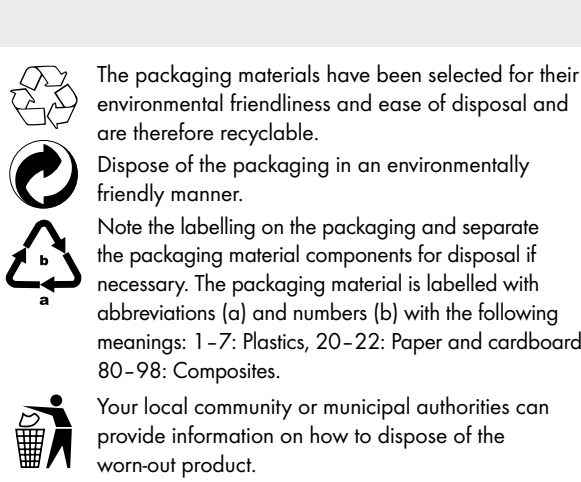Ordering replacement parts

You can order replacement brush heads for this product on the internet at **www.kompernass.com**.

- Scan the QR code with your smartphone/tablet. You can use this QR code to go directly to our website to view and order the available spare parts. If you have problems with your online order, you can contact our service centre by phone or e-mail. Please always quote the article number (e.g. 123456_7890) when ordering. You can find it on the title page of these operating instructions.

Kompernass Handels GmbH warranty

Dear Customer, This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred. If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

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- 15 -

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure
To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 373276_2104 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

- You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com . This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 373276_2104.

Service
GB Service Great Britain Tel.: 0800 044 7657 E-Mail: kompernass@lidl.co.uk
MT Service Malta Tel.: 80062230 E-Mail: kompernass@lidl.com.mt IAN 373276_2104

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

KOMPERNASS HANDELS GMBH
BURGSTRASSE 21 - 44867 - BOCHUM GERMANY
www.kompernass.com

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■ GB/MT

Introduction

Congratulations on the purchase of your new appliance. You have selected a high-quality product. The operating instructions are part of this product. They contain important information on safety, usage and disposal. Before using the product, familiarise yourself with all handling and safety guidelines. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner (s).

Intended use

This appliance is intended exclusively for personal dental care. This appliance is intended solely for use in private homes. The appliance is not intended for use in commercial, industrial or medical environments. It is not suitable for animal care.

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■ Children must not use the appliance as a playing.
■ The appliance must not come into contact with hot surfaces. The appliance could be irreparably damaged!
■ Do not open the appliance (except the battery compartment) yourself. Risk of injury if appliance is not properly assembled!
Information on using batteries
Please observe the following when handling batteries:
■ Do not throw batteries into a fire. Do not recharge batteries. There is a risk of explosion and injury!
■ Never open batteries and never solder or weld batteries. There is a risk of explosion and injury!
■ Check the condition of the batteries at regular intervals. Leaking batteries can cause damage to the appliance.
■ If you do not intend to use the appliance for an extended period, remove the batteries.
■ If the batteries have leaked, use protective gloves. Clean the battery compartment and the battery contacts with a dry cloth.
■ Keep batteries out of the reach of children. Children might put batteries into their mouth and swallow them. If a battery is swallowed, seek medical assistance immediately.
■ Batteries must be inserted with the correct polarity.
■ Always use the stated battery type.
■ Do not use different types of batteries together or mix new batteries with used batteries.
■ Do not short-circuit the connection terminals.
■ Remove depleted batteries from the appliance and dispose of them safely.
Operation
NOTE
► Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.
■ Remove the battery compartment cover ❹ from the hand element ❸ with a slight tilting movement and insert the batteries according to the polarity diagram in the battery compartment. Close the cover of the battery compartment ❹ again.

- 19 -

■ Children must not use the appliance as a playing.
■ The appliance must not come into contact with hot surfaces. The appliance could be irreparably damaged!
■ Do not open the appliance (except the battery compartment) yourself. Risk of injury if appliance is not properly assembled!
Information on using batteries
Please observe the following when handling batteries:
■ Do not throw batteries into a fire. Do not recharge batteries. There is a risk of explosion and injury!
■ Never open batteries and never solder or weld batteries. There is a risk of explosion and injury!
■ Check the condition of the batteries at regular intervals. Leaking batteries can cause damage to the appliance.
■ If you do not intend to use the appliance for an extended period, remove the batteries.
■ If the batteries have leaked, use protective gloves. Clean the battery compartment and the battery contacts with a dry cloth.
■ Keep batteries out of the reach of children. Children might put batteries into their mouth and swallow them. If a battery is swallowed, seek medical assistance immediately.
■ Batteries must be inserted with the correct polarity.
■ Always use the stated battery type.
■ Do not use different types of batteries together or mix new batteries with used batteries.
■ Do not short-circuit the connection terminals.
■ Remove depleted batteries from the appliance and dispose of them safely.
Operation
NOTE
► Due to the unaccustomed