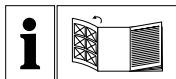
MASAJEADOR ANTICELULÍTICO / DISPOSITIVO PER MASSAGGI ANTICELLULITE SMG 2.4 A1

ES	IT MT
MASAJEADOR ANTICELULÍTICO	DISPOSITIVO PER
Instrucciones de uso	MASSAGGI ANTICELLULITE
	Istruzioni per l'uso
PT	GB MT
MASSAJADOR DE CELULITE	CELLULITE MASSAGER
Manual de instruções	Operating instructions
DE AT CH	
CELLULITE-MASSAGEGERÄT	
Bedienungsanleitung	

IAN 385372_2107

ES

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www.kompernass.com

Contents

Introduction	62
Intended use	62
Package contents	63
Disposal of the packaging	64
Appliance description	65
Technical data	66
Safety information	68
Use	73
Cleaning	74
Storage	75
Disposal of the appliance	75
Mains adapter disposal	75
Kompernass Handels GmbH warranty	76
Service	79
Importer	79

Introduction

Congratulations on the purchase of your new appliance.
You have chosen a high-quality product.

The operating instructions are an integral part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

Intended use

This appliance is designed for massaging the skin on the buttocks, hips, upper thighs and upper arms. The massaging movements of the appliance promote the blood circulation in the corresponding skin areas.

This appliance is intended solely for personal, private use, it is not to be used for medical or commercial applications and is not a substitute for medical treatment. The appliance is designed for use in dry indoor rooms. It is not intended for use in the medical/therapeutic or commercial fields.

 WARNING**Danger if not used as intended!**

Failure to use the appliance for its intended purpose and/or any other type of use can be hazardous.

- ▶ Use the appliance only for its intended purpose.
- ▶ Observe the procedures described in these operating instructions.

Package contents

The appliance is supplied with the following components as standard:

- Cellulite massager
 - Mains adapter
 - Operating instructions
- 1) Remove all parts of the appliance and the operating instructions from the box.
 - 2) Remove all packaging material.

NOTE

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).

Disposal of the packaging

The packaging protects the appliance from damage during transport. The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable.

Returning the packaging into the material cycle conserves raw materials and reduces the amount of waste generated. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1–7: Plastics, 20–22: Paper and cardboard, 80–98: Composites.

NOTE

- If possible, keep the original packaging during the warranty period so that the appliance can be properly packed for returning in the event of a warranty claim.

Appliance description

Figure A:

- ❶ Adjustable handle
- ❷ Switch: **0** = OFF, **1** = Slow rotation, **2** = Fast rotation
- ❸ Socket
- ❹ 6 removable massage rolls

Figure B:

- ❺ Plug
- ❻ Mains adapter

Technical data

Mains adapter	
Manufacturer	Shenzhen Sunlight Electronic Technology Co., Ltd., China EU representative: EUREP GmbH, Unterlettenweg 1a, 85051 Ingolstadt, Germany, District Court Ingolstadt HRB 9948
Model	D03DE1200120A
Input voltage	100–240 V ~ (AC)
Input AC frequency	50/60 Hz
Output voltage	12.0 V $\overline{=}$ (DC)
Output current	1.2 A
Output power	14.4 W
Average operating efficiency	84.9%
Low load efficiency (10%)	80.5%
No-load power consumption	0.05 W
Input current	0.7 A
Efficiency class 6	Ⓥ

Mains adapter

Polarity	
Safety transformer, short-circuit proof	
Switched-mode power supply	
Protection class	II /  (double insulation)
Nominal ambient temperature (ta)	40 °C
Protection type	IP20: Protection against ingress of solid objects with a diameter of more than 12.5 mm

Appliance

Input voltage	12.0 V 
Input current	1.2 A
Polarity	

Safety information

RISK OF ELECTRIC SHOCK!

- ▶ Operate the appliance only with the supplied mains adapter.
- ▶ The appliance should only be connected to correctly installed and earthed mains power sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- ▶ Ensure that the power cable does not become wet or moist during use. Route the cable in such a way that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs to the appliance. All repairs must be performed by our Customer Service Unit or by a qualified specialist technician.
- ▶ Always unplug the appliance when it is not in use and before you clean it.

** RISK OF ELECTRIC SHOCK!**

- ▶ Never lift, carry or drag the appliance by the power cable.
- ▶ To prevent accidents, have defective plugs, power cables and/or mains adapters replaced immediately by an authorised specialist, our Customer Service or a similarly qualified person.
- ▶ Use the appliance indoors only; never use it outside.
- ▶ Do not use the appliance while taking a bath or shower. Do not store the appliance anywhere where it could fall into the bath tub or the wash basin.
- ▶ Never try to retrieve an appliance that has fallen into water. Remove the mains adapter from the power socket first.
- ▶ Do not use the appliance if the power cable is damaged, the mains adapter is defective, the appliance is not working properly, if it has fallen into water or is otherwise damaged.

** RISK OF ELECTRIC SHOCK!**

- ▶ Keep the power cable away from the rotating massage rollers.

Do not use the mains adapter if it has a damaged mains plug or cable.

** WARNING! RISK OF INJURY!**

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- ▶ Do not allow children to use the appliance as a toy.

** WARNING! RISK OF INJURY!**

- ▶ People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- ▶ Do not use the appliance if you suffer from one or more of the following illnesses or complaints: circulatory problems, varicose veins, open wounds, bruising, skin tears, inflammation of the veins (phlebitis) or thrombosis.
- ▶ Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.
- ▶ Do not use the appliance under blankets or pillows.
- ▶ Do not use the appliance after taking substances that could reduce your capacity to react quickly (e.g. painkillers or alcohol).

** WARNING! RISK OF INJURY!**

- ▶ Stop treatment with the appliance immediately if you experience any discomfort or pain.
- ▶ Uninterrupted application of the appliance should not last any longer than 15 minutes. Afterwards, switch off the appliance for at least half an hour and let it cool down.
- ▶ Do not put your fingers or any other objects between the rotating massage rollers.
- ▶ Use the appliance only for its intended purpose.

CAUTION! PROPERTY DAMAGE!

Use the mains adapter/the appliance indoors only.

- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surface.

NOTE

- ▶ No user action is required to switch the product between 50 and 60 Hz. The product adapts automatically to either 50 or 60 Hz.

Use

WARNING! RISK OF INJURY!

- ▶ If you have long hair, tie it back before using the appliance on your upper arms. Otherwise, your hair could get caught in the appliance.
- ▶ Please beware of your long hair winded into the massage rolls **4** while using the machine.
- ◆ Insert the plug **5** on the mains adapter **6** into the socket **3** on the massager and then insert the mains plug **6** into a wall outlet.
- ◆ Slide your hand under the adjustable handle **1**. If this is too wide or too narrow, you can adjust it to fit your hand using the Velcro fastener.
- ◆ Press the massager gently against the part of the body that you want to massage.
- ◆ Switch on the appliance using the switch **2**. Select position **1** for a slow rotation of the rollers and position **2** for a fast rotation.

NOTE

- ▶ If you find the application unpleasant, try reducing the speed, adjust the pressure applied, change to a different part of the body or stop the massage completely.
- ◆ To switch the appliance off, push the switch **2** back to the **0** position.

Cleaning

** RISK OF ELECTRIC SHOCK**

- ▶ Before cleaning the appliance, always pull out the mains adapter **6** from the mains power socket.

CAUTION! PROPERTY DAMAGE!

- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surfaces.
- For easier cleaning, pull the massage rollers **4** off of the appliance. Clean them in warm water with a little mild soap. Rinse off the soap residues thoroughly and allow the massage rollers **4** to dry in the air. When they are dry, push them back onto the drive shaft of the massager.
- Wipe the appliance with a damp cloth. If necessary, add a little detergent to the cloth and wipe the appliance off afterwards with a cloth dampened with clear water.
- Wipe the mains adapter **6** and mains cable with a damp cloth. Make sure the mains adapter **6** has completely dried off before the next use.

Storage

- Clean the appliance as described in the section "Cleaning".
- Store the appliance in a clean and dust-free location.

Disposal of the appliance

Never dispose of the appliance in your normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance through an approved disposal facility or your municipal waste facility. Please observe the currently applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Mains adapter disposal

Do not dispose of the mains adapter in your normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the mains adapter via an approved disposal centre or at your community waste facility. Please observe the currently applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Kompernass Handels GmbH

warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 385372_2107 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 385372_2107.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

MT Service Malta

Tel.: 80062230

E-Mail: kompernass@lidl.com.mt

IAN 385372_2107

Importer

Please note that the following address is not the service address.
Please use the service address provided in the operating instructions.

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