PDF ONLINE
www.lidl-service.com

WATER JET FLOSSER SMD 24 B2

(GB) (IE) (NI) (CY)

WATER JET FLOSSER

Operating instructions and safety instructions

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
IAN 444782_2307

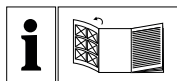GB IE NI CY

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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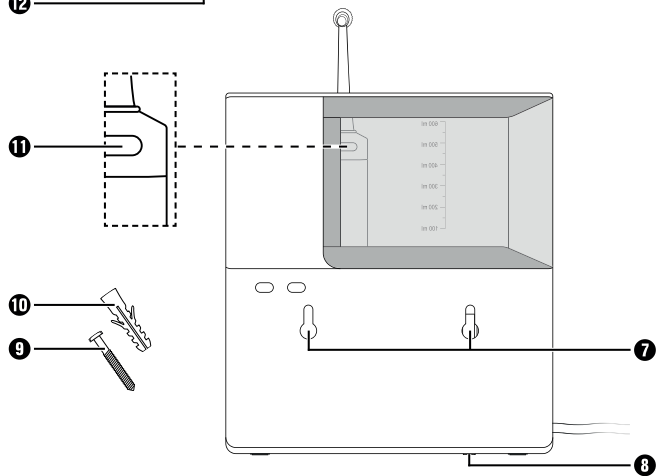
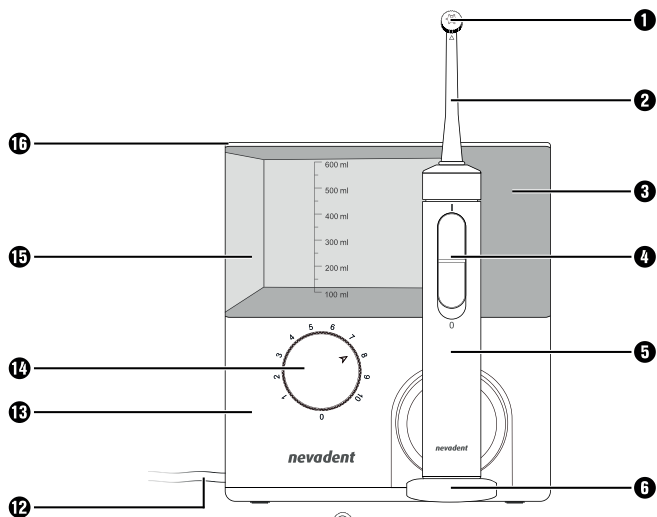
Contents

Introduction. 2

Information about these operating instructions 2
Intended use. 2

Warnings and symbols used 3

Safety information 4

Package contents and transport inspection 7

Description of the appliance 8

Before first use 9

Requirements for the set-up location. 9
Wall mounting 9

Handling and operation. 11

Attaching/removing nozzle attachments . 11
Adjusting the nozzle control. 12
Using the oral irrigator. 13
After use. 15

Cleaning and care. 16

Storage 17

Troubleshooting 18

Ordering replacement parts . . . 20

Disposal. 21

Disposal of the appliance. 21
Disposal of the packaging 21

Appendix. 22

Technical data 22
Kompernass Handels GmbH warranty . . 23
Service 25
Importer 25

Introduction

Information about these operating instructions

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

Intended use

This appliance is intended exclusively for dental care and mouth hygiene in humans. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. This appliance is intended solely for use in private households. It is not intended for any other purpose nor for use beyond the scope described.

Claims of any kind for damage resulting from misuse, incompetent repairs, unauthorised modification or the use of unauthorised replacement parts will not be accepted. The operator bears sole liability.

Warnings and symbols used

The following warnings and symbols are used in these operating instructions, on the packaging and on the appliance (where applicable):

	DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminently hazardous situation that will result in death or serious injury if not avoided.
	WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation that could result in death or serious injury if not avoided.
	CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation that could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation that could result in property damage if not avoided.
	A note identifies additional information that facilitates the use of the appliance.
	AC current/voltage

Safety information

⚠ DANGER! ELECTRIC SHOCK!

- ▶ Connect the appliance only to a properly installed and easily accessible mains power socket supplying a mains power voltage of 100–240 V ~, 50/60 Hz.
- ▶ Ensure that the power cable is not damaged. Keep it out of hot areas and route it in such a way that it cannot be damaged.
- ▶ Before use, check the appliance for visible external damage. Never use the appliance if the cable, plug or housing are damaged.
- ▶ To avoid potential risks, if the appliance power cable is damaged, it must be replaced by the manufacturer, its customer service or by a qualified technician.
- ▶ Have all repairs carried out by a specialist workshop. Under no circumstances should you open the appliance yourself. Repairs that are not carried out by a specialist workshop could lead to physical injuries.
- ▶ Do not use the appliance in the vicinity of a bathtub/shower or a sink filled with water. The proximity of moisture presents a danger even when the appliance is switched off.

- ▶ Never touch the appliance if it has fallen into water. Always disconnect the power plug from the power socket first!
- ▶ Never immerse the connected appliance in water or other liquids.
- ▶ Never touch the power cable or the appliance with wet or damp hands.
- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way.
- ▶ Do not use the power cable with an extension cable; connect the power cable directly to a power socket.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not change any accessories while the appliance is switched on.
- ▶ Always store the appliance indoors. To prevent accidents, keep the appliance in a dry location when not in use.
- ▶ Oral irrigators may be used by children over the age of 8 and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised. Do not allow children to play with the appliance.

- ▶ Switch off the appliance before you put it down or change accessories.
- ▶ Lay the power cable so that nobody can step on it or trip over it.
- ▶ Always pull out the plug before cleaning the appliance. Do not pull on the power cable. Always grip the power plug to disconnect the appliance from the power supply.

! ATTENTION! PROPERTY DAMAGE!

- ▶ Use only accessories recommended by the manufacturer.
- ▶ Do not make any unauthorised modifications or alterations to the appliance.
- ▶ Never put down the appliance next to radiators, ovens or other heated appliances or surfaces.
- ▶ The appliance may only be used with cold or lukewarm tap water up to a maximum of 40°C . If necessary, add a few drops of mouthwash. The water tank should never be completely filled with additives such as mouthwash or mouth rinse. These can damage the appliance.

Package contents and transport inspection

- ◆ Remove all parts of the appliance and the operating instructions from the carton.
- ◆ Remove all packaging materials and any films and labels.

DANGER!

- Do not allow children to play with packaging materials. There is a risk of suffocation!

The package contents include the following components (see fold-out page for illustrations):

- Water jet flosser (main unit with hand unit and water tank)
- 4 nozzle attachments
- 2 wall plugs (ø 5 mm)
- 2 screws (3 mm)
- Operating instructions (not shown)

Note

- Check the package for completeness and for signs of visible damage.
- If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the service hotline (see section **Service**).

Description of the appliance

(See fold-out page for illustrations)

- ❶ Nozzle control
- ❷ Nozzle attachment
- ❸ Storage compartment for nozzle attachments
- ❹ On/Off switch for the water jet I/O
- ❺ Hand unit
- ❻ Cradle for hand unit
- ❼ Hangers
- ❽ Locking pin
- ❾ Screw
- ❿ Wall plug
- ⓫ Release button for nozzle attachment
- ⓬ Power cable (with mains plug)
- ⓭ Main unit
- ⓮ Rotary switch (On/Off switch and water pressure control)
- ⓯ Water tank
- ⓰ Lid

Before first use

Requirements for the set-up location

To ensure safe and trouble-free operation of the appliance, the set-up location must fulfil the following requirements:

- The appliance may only be operated free-standing if used in dry rooms. The set-up location must be a horizontal, level and stable surface.
- In damp rooms, especially in the vicinity of water, the appliance must be firmly fixed to the wall to protect it from falling into water.
- The mains power socket must be easily accessible so that the power cable can be quickly disconnected in an emergency.

Wall mounting

Note

- Before mounting to the wall, make sure that there are no electric cables, wires or pipes running behind the wall surface where you plan to hang the appliance.
1. Mark two drill holes on the wall 7.2 cm apart and drill the marked holes with a 5 mm bit. Ensure that there is a power socket in the immediate vicinity.
 2. Push the two supplied wall plugs **10** into the drilled holes and screw the two supplied screws **9** into the wall plugs **10**. Take care that the screws **9** still protrude about 3 mm out of the holes.

3. On the back side of the appliance, the right hanger **7** has a locking pin **8** to prevent the appliance from coming free of the screws **9**. Pull the locking pin **8** down slightly out of the hanger **7** (see Fig. 1).

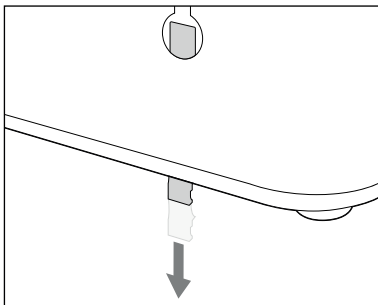
Fig. 1

4. Push the appliance onto the screws **9** so that the screws **9** enter into the hangers **7** on the back side of the appliance. Then press the appliance downward so that the screws **9** slide into the rails of the hangers **7** and the appliance rests firmly on the screws **9**.
5. Slide the locking pin **8** upwards until it audibly engages with the hanger **7**.
6. Plug the mains plug of the power cable **12** into a socket.

Handling and operation

Attaching/removing nozzle attachments

i Note

- ▶ If used daily, replace the nozzle attachments **2** after no more than six months.
- ▶ Use only original nozzle attachments **2**. Information on ordering replacement nozzle attachments **2** can be found in the section **Ordering replacement parts**.
- ◆ Make certain that the appliance is switched off before attaching or removing the nozzle attachment **2**.
- ◆ Fit a nozzle attachment **2** onto the hand unit **5**. Press the nozzle attachment **2** into the opening until the nozzle attachment **2** audibly engages in the hand unit **5**. Ensure that the guide tab on the lower end of the nozzle attachment **2** is aligned with the corresponding groove in the opening on the hand unit **5** (see Fig. 2).
- ◆ To remove the nozzle attachment **2**, press the release button **11** on the hand unit **5** and then pull the nozzle attachment **2** out of the hand unit **5** (see Fig. 3).

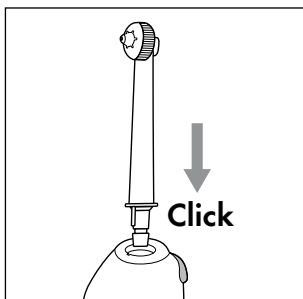
Fig. 2

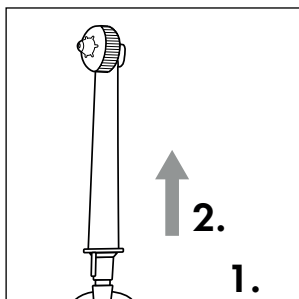
Fig. 3

Adjusting the nozzle control

You can use the nozzle control **1** on the head of the nozzle attachment **2** to switch between two different nozzle functions:

- ◆ Turn the nozzle control **1** anti-clockwise with your fingers to obtain a single jet for targeted removal of food remnants from hard-to-reach areas (see Fig. 4 and 5).

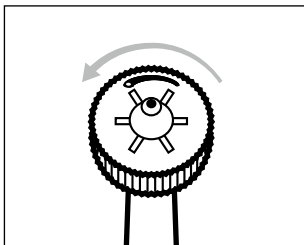
Fig. 4

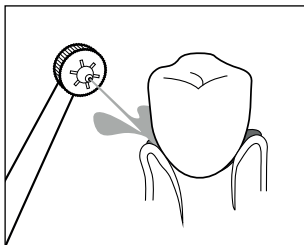
Fig. 5

- ◆ Turn the nozzle control **1** clockwise with your fingers to obtain multiple jets for massaging the gums (see Fig. 6 and 7).

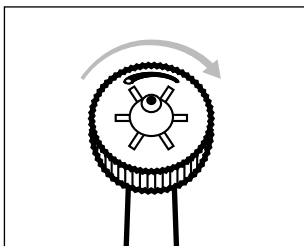
Fig. 6

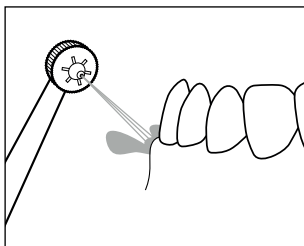
Fig. 7

Using the oral irrigator

CAUTION!

- ▶ Fill the water tank **15** with pure tap water. If necessary, add a few drops of mouthwash. The water tank **15** should never be completely filled with additives such as mouthwash or mouth rinse.
- ▶ Use only cold or lukewarm water with a maximum temperature of 40°C. Do not fill the water tank **15** with hot water.

Note

- ▶ Use of the appliance is not a replacement for daily tooth cleaning with a toothbrush and toothpaste. Use the appliance only as a supplement to daily tooth cleaning with toothbrush and toothpaste.
- ▶ Only use the appliance for a maximum of 2 minutes or only long enough to use up one full tank of water (600 ml). One full tank of water (600 ml) is sufficient to operate the appliance for about 150 seconds at level **1** or about 110 seconds at level **10**.
- ▶ The unfamiliar stimulation may cause some gum bleeding at first, which is not unusual. If this continues for more than two weeks, you should consult your dentist.
- ▶ Consult your dentist before using the appliance if you suffer from severe periodontitis, have injuries or abscesses in the mucous membrane of the oral cavity or underwent dental surgery within the last two months.
- ◆ Remove the water tank **15** from the main unit **13** and open the lid **16**. Fill the water tank **15** with cold or lukewarm water. Take care not to fill the water tank **15** past the maximum fill mark of **600 ml**. Then close the water tank **15** again with the lid **16** and place the water tank **15** in the main unit **13**.

- ◆ Fit a nozzle attachment **2** onto the hand unit **5**. Press the nozzle attachment **2** into the opening until the nozzle attachment **2** audibly engages in the hand unit **5**. Ensure that the guide tab on the lower end of the nozzle attachment **2** is aligned with the corresponding groove in the opening on the hand unit **5**.
- ◆ Switch on the main unit **13** by setting the rotary switch **14** to the desired water pressure level **1 - 10** (**1** = gentle, **10** = strong).
- ◆ Open your mouth slightly and insert the nozzle attachment **2**. Bend over a washbasin to allow the water to run into the basin.

i Note

- ▶ Begin with a lower level to allow your gums to become accustomed to the water jet. Increase the strength of the water jet in stages. If the pressure is uncomfortable, select a lower level again.
- ▶ To avoid splashes, do not switch on the hand unit **5** until it is in your mouth.
- ◆ Switch on the hand unit **5** by moving the On/Off switch for the water jet **I/O 4** on the hand unit **5** to **I**.
- ◆ Direct the water jet at the gaps between your teeth and the edges of your teeth to clean them. Position the nozzle attachment **2** so that the water jet can rinse the food remnants away from the gums (see Fig. 8). Move the nozzle attachment **2** slowly along the edge of the teeth from gap to gap (see Fig. 9). Clean all the front and back sides of your teeth in this way.

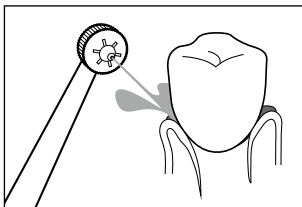
Fig. 8

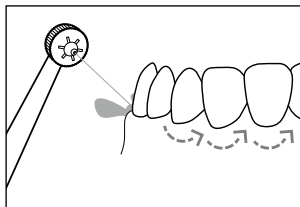
Fig. 9

- ◆ You can use the image below as a guide for the order of cleaning (see Fig. 10). The numbers in the image indicate the order in which to clean the different areas:

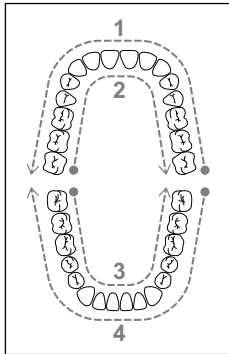
Fig. 10

1. Front of upper teeth
2. Back of upper teeth
3. Back of lower teeth
4. Front of lower teeth

- ◆ Start at the front side of your upper molars and move the nozzle attachment **2** along the edges of the teeth from one side to the other. Repeat this process for the back of the upper teeth and for the back and front of the lower teeth.

After use

i Note

- Empty the water tank **15** after every use to prevent the growth of germs and bacteria.
1. Switch the On/Off switch for the water jet **I/O 4** on the hand unit **5** to **0** before removing the nozzle attachment **2** from your mouth.
 2. Set the rotary switch **14** on the main unit **13** to **0**.

3. Remove the water tank **15** from the main unit **13** and open the lid. Pour the remaining water out of the water tank **15** into the washbasin.
4. Close the lid **16** and place the water tank **15** back in the main unit **13**.
5. Set the rotary switch **14** on the main unit **13** to **10** and hold the hand unit **5** over the washbasin.
6. Switch the On/Off switch for the water jet I/O **4** on the hand unit **5** to **I** to clear any remaining water from the line and rinse out the nozzle attachment **2**.
7. Switch the On/Off switch for the water jet I/O jet **4** on the hand unit **5** to **0**.
8. Set the rotary switch **14** on the main unit **13** to **0**.
9. Press the release button **11** on the hand unit **5** and then pull the nozzle attachment **2** out of the hand unit **5**.
10. Rinse the nozzle attachment **2** under running water.
11. Dry the nozzle attachment **2** and place it in one of the spots in the storage compartment **3**.
12. Place the hand unit **5** in the cradle **6** on the main unit **13**.

Cleaning and care

DANGER!

- ▶ Always pull out the plug before cleaning the appliance. There is a risk of electric shock!
- ▶ The appliance must be completely dry before connecting it to the socket again.

Never submerge the main unit **13** or the hand unit **5** in water or other liquids!

CAUTION!

- ▶ Do not use any aggressive cleaning agents or solvents. These could damage the plastic surfaces of the appliance.
- ◆ Remove the nozzle attachment **2** from the hand unit **5** and rinse it off under running water.
- ◆ Clean the main unit **13** and the hand unit **5** with a soft, slightly damp cloth. For stubborn soiling add a little mild detergent to the cloth. Afterwards, wipe off with a cloth moistened only with pure water. Then dry off the main unit **13** and the hand unit **5** with a dry cloth.
- ◆ Remove the water tank **15** from the main unit **13** and open the lid **16**. Rinse the water tank **15** under running water. After cleaning, pour the remaining water out of the water tank **15** into the washbasin.
- ◆ If you do not intend to use the appliance for a longer period, remove the lid **16** from the water tank **15**. Allow the water tank **15** to dry completely in air before placing it back in the main unit **13**.

Storage

- ◆ Disconnect the appliance from the power supply by pulling the plug of the switched-off appliance from the socket.
- ◆ Remove any water from the appliance and allow the appliance to dry completely (see section **Cleaning and care**).
- ◆ Store the appliance in a dry and dust-free location out of direct sunlight.

Troubleshooting

The following table will help you to identify and rectify minor malfunctions. If the solutions presented below do not resolve the problem, please contact Customer Service (see section **Service**).

Fault	Cause	Remedy
No water jet comes out of the nozzle attachment 2 .	The appliance is not switched on.	Make sure that the appliance is properly connected to the power supply. Set the rotary switch 14 to the desired water pressure level 1 – 10 . Switch the On/Off switch for the water jet I/O 4 on the hand unit 5 to I.
	The mains plug is not plugged in.	Insert the plug into a mains power socket.
	The socket is not supplying any power.	Check the house mains fuse.
	The nozzle attachment 2 is defective.	Replace the nozzle attachment 2 (see section Attaching/removing nozzle attachments)
	The water tank 15 is empty.	Fill the water tank 15 with tap water (see section Using the oral irrigator).
	The appliance is defective.	Switch off the appliance and pull out the mains plug. Contact Customer Service (see section Service).

Fault	Cause	Remedy
The water pressure cannot be set correctly or is too weak.	The water tank 15 is not seated properly in the main unit 13 .	Place the water tank 15 in the main unit 13 and press the water tank 15 down until it sits firmly in the main unit 13 .
The water tank 15 has a leak.	The water tank 15 is not seated properly in the main unit 13 .	Place the water tank 15 in the main unit 13 and press the water tank 15 down until it sits firmly in the main unit 13 .
	The valve of the water tank 15 is defective.	Switch off the appliance and pull out the mains plug. Contact Customer Service (see section Service).
The nozzle attachment 2 has a leak.	The nozzle attachment 2 is not properly attached to the hand unit 5 .	Attach the nozzle attachment 2 properly to the hand unit 5 (see section Attaching/removing nozzle attachments).
	The nozzle attachment 2 is defective.	Replace the nozzle attachment 2 (see section Attaching/removing nozzle attachments).

Ordering replacement parts

You can order replacement parts for this product on the Internet at **www.kompernass.com**.

Scan this QR code with your smartphone or tablet.
You can use this QR code to go directly to our website **www.kompernass.com** to view and order the available replacement parts.

Note

- ▶ If you have problems with your online order, you can contact our service centre by phone or e-mail.
- ▶ Always quote the article number (IAN) **444782_2307** in your order.
- ▶ Please note that online ordering of replacement parts is not possible for all countries.

Disposal

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

1–7: plastics, 20–22: paper and cardboard, 80–98: composites.

Appendix

Technical data

Power supply	100–240 V ~ (AC), 50/60 Hz
Power consumption	24 W
Protection class	II/□ (double insulation)
Protection type	IPX4 (protected against water splashes from any direction)

Note

- No user action is required to switch the product between 50 and 60 Hz.
The product switches automatically to either 50 or 60 Hz.

This mark confirms that the product complies with the product safety requirements applicable in the UK.

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 444782_2307 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 444782_2307.

Service

GB

Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE

Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

CY

Service Cyprus

Tel.: 8009 4241

E-Mail: kompernass@lidl.com.cy

IAN 444782_2307

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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Lidl Great Britain Ltd