WATER JET FLOSSER / MUNDDUSCHE HYDROPULSEUR CBMD 3.7 A1

(GB) (IE) (NI) (MT)

WATER JET FLOSSER

Operating instructions

(FR) (BE) (CH)

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Mode d'emploi

(IT) (CH) (MT)

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(PT)

HIDROPULSOR PARA HIGIENE ORAL

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(DE) (AT) (BE) (CH)

MUNDDUSCHE

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(NL) (BE)

MONDDOUCHE

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
(ES)

IRRIGADOR BUCAL

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IAN 500431_2504

(GB) (IE) (NI) (DE) (AT) (BE)
(CH) (FR) (NL) (ES) (IT) (MT) (PT)

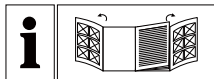GB IE NI MT

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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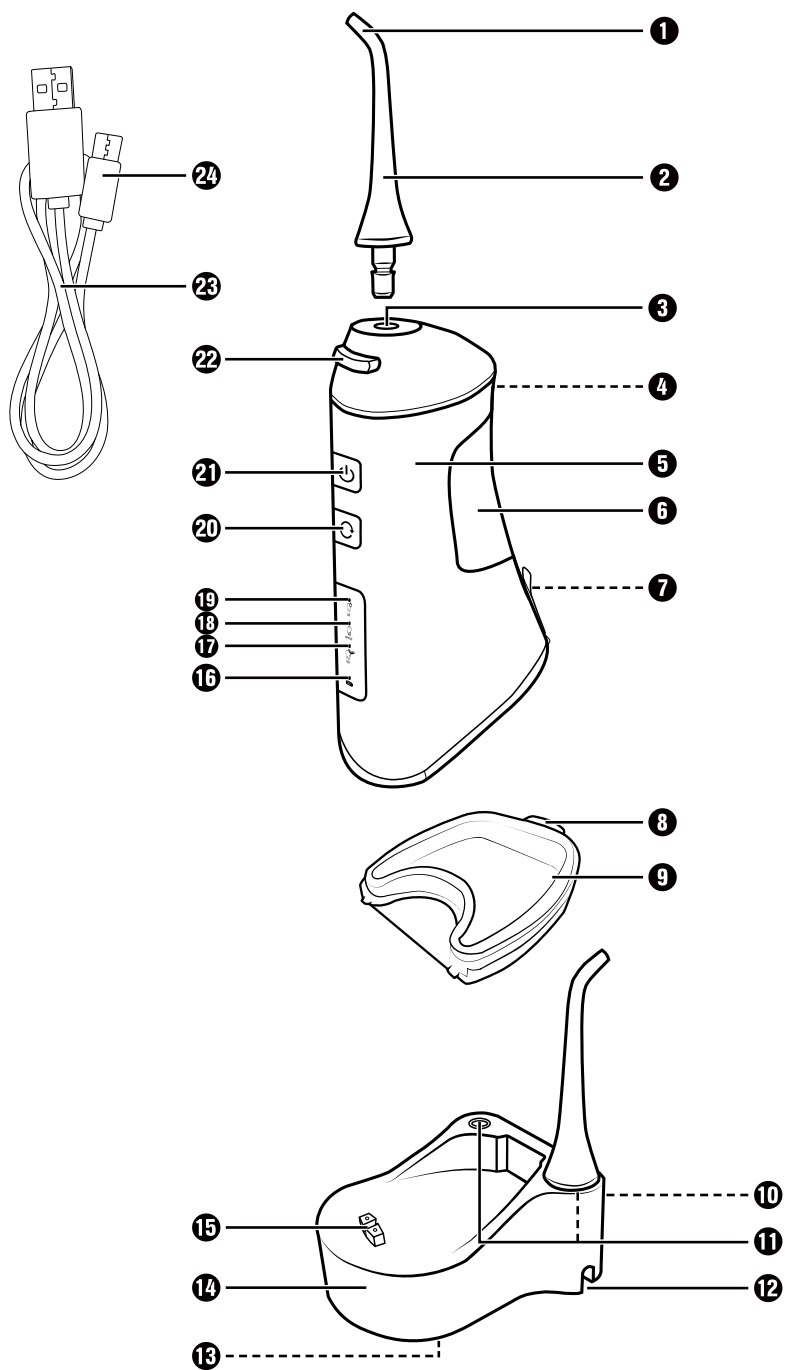
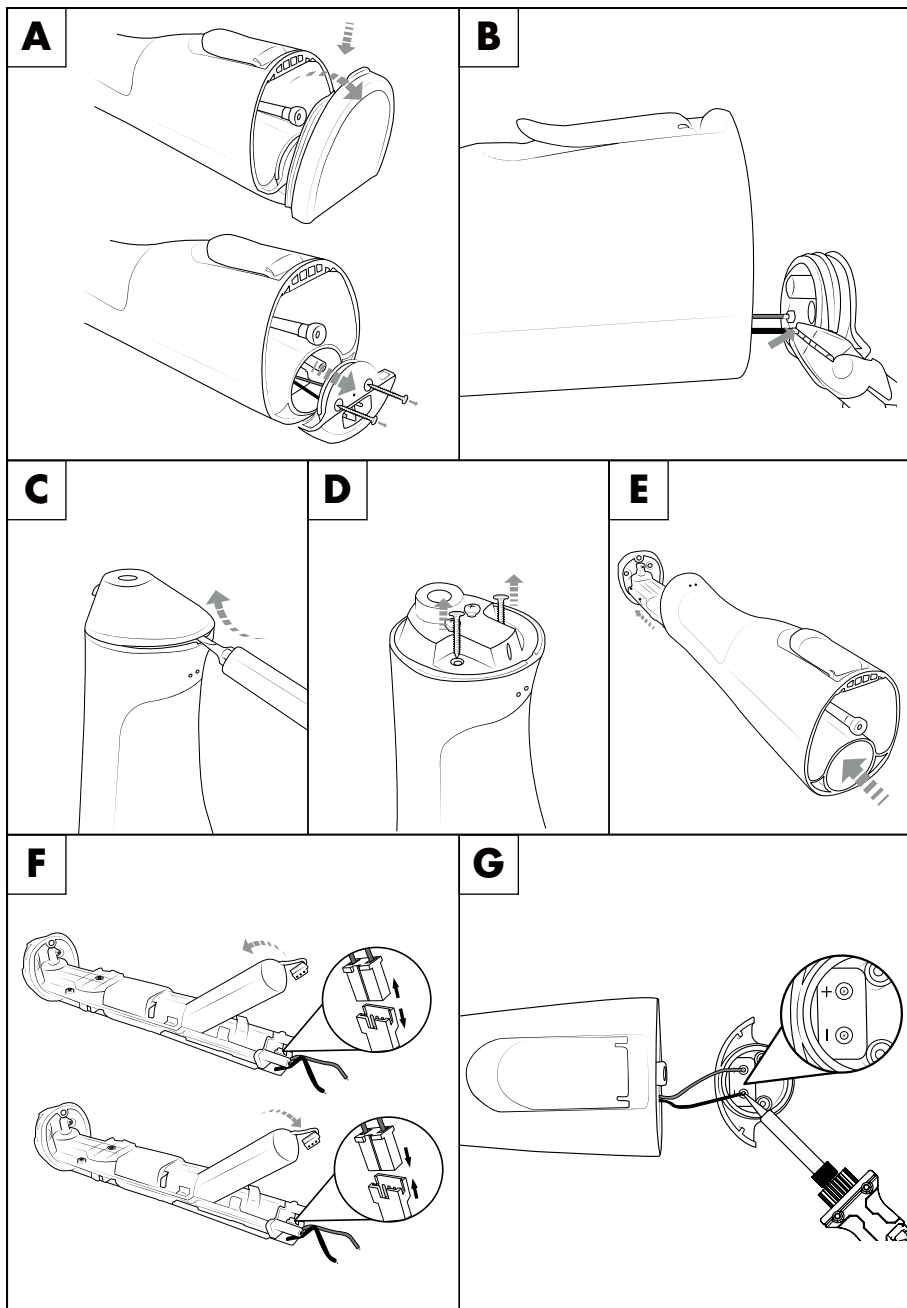
Contents

1. Information about these operating instructions	2
2. Intended use	2
3. Warnings and symbols used	2
4. Safety information	3
5. Technical data	5
6. Package contents	5
7. Appliance description	6
8. Setting up the appliance / wall mounting	6
9. Charging the appliance	7
10. Filling the water tank	7
11. Fitting / removing attachments	7
12. Before first use	8
13. Power levels	8
14. Operation	8
15. After use	10
16. Cleaning and maintenance	10
17. Storage	10
18. Replacing the battery	11
19. Disposal	12
19.1. Disposal of the appliance and packaging	12
19.2. Disposal of batteries	12
20. Kompernass Handels GmbH warranty	12
20.1. Service	14
20.2. Importer	14
21. Ordering replacement parts	14

1. Information about these operating instructions

Congratulations on the purchase of your new appliance. You have selected a high-quality product.

The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass on all documents to any future owner.

2. Intended use

This appliance is intended exclusively for personal dental care at home. It is intended solely for use in private households. This appliance is not intended for use in industrial, commercial or medical environments. It is not suitable for animal care. It is not intended for any other purpose nor for use beyond the scope described.

3. Warnings and symbols used

The following warnings and symbols are used in these operating instructions, on the packaging and on the device (where applicable):

	DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminently hazardous situation that will result in death or serious injury if not avoided.
	WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation that could result in death or serious injury if not avoided.

	CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation that could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation that could result in property damage if not avoided.
	Note: A note identifies additional information that facilitates the use of the appliance.
	DC current/voltage
	Removable power supply
	Do not dispose of electrical appliances in the household waste!
	Dispose of the packaging in an environmentally-responsible manner.
	The packaging is made of recyclable materials. Observe the labelling on the packaging material when sorting for disposal: The material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

  ES/PT	The packaging contains paper and/or cardboard components. The packaging contains plastic and/or metal components.
UK CA	This mark confirms that the product complies with the product safety requirements applicable in the UK.

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BAC OR TBI

FR: The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

4. Safety information

WARNING! RISK OF INJURY!

- Check the appliance for visible external damage before use. Never use the appliance if the charging cable, plug or housing are damaged.
- In the event of malfunctions or obvious damage, please contact our Customer Service department.
- Have all repairs carried out by a specialist workshop. Under no circumstances should you open the appliance yourself. Repairs that are not carried out by a specialist workshop could lead to physical injuries.
- Do not use the appliance in the vicinity of a bathtub/shower or a sink filled with water. The proximity of moisture presents a danger even when the appliance is switched off.
- To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- Ensure that no foreign objects enter the appliance.
-  Do not submerge the handpiece or charging station in water or other liquids.

- This appliance may be used by children from the age of 8 and by persons with reduced physical, sensory or mental abilities or lack of experience and/or knowledge if they are supervised.
- Do not allow children to play with the appliance. Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- Always store the appliance indoors to avoid accidents. After use, store the appliance in a dry place away from direct sunlight.
- Under extreme conditions, leakage can occur from the battery cells and electrolyte solution can leak out. If the liquid comes into contact with skin or eyes, immediately rinse the affected area with plenty of clean water and seek medical attention.
- Switch off the appliance before you put it down or change accessories.

⚠ ATTENTION! PROPERTY DAMAGE!

- Charge the appliance only with the supplied charging station.
- Use only accessories recommended by the manufacturer.
- Never put down the appliance next to or on radiators, ovens or other heated appliances or surfaces.
- Ensure that the charging cable is not stretched tightly or kinked.
- The appliance may only be used with cold or lukewarm tap water up to a maximum of 40°C. If necessary, add a few drops of mouthwash. The water tank should never be completely filled with additives such as mouthwash or mouth rinse. These can damage the appliance.
- This device contains a rechargeable battery that can only be replaced by trained personnel.

- The product is equipped or supplied with a lithium-ion battery. Do not throw the rechargeable battery into a fire and do not subject it to high temperatures. There is a risk of an explosion!
-  Do not charge or use the appliance outdoors.
- Never store the device when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the device is stored for a longer period of time, the charge level of the battery must be checked regularly. The optimum charge level is between 50% and 80%.

5. Technical data

Charging station	
Input voltage/ current	5 V  , 2 A (via the USB-C port)
Output voltage/ current	5 V  , 2 A
Protection class	III /  (protection through low voltage)
Battery	
Capacity	1200 mAh
Battery (lithium ion)	3.7 V  , 4.44 Wh
Handpiece	
Input voltage/ current	5 V  , 2 A
Protection class	III /  (protection through low voltage)
Protection class	IPX7 (protection against temporary submersion)
Water tank volume	210 ml
Operating time with full battery charge	Approx. 20 minutes

6. Package contents

The appliance is supplied with the following components as standard:

- Handpiece with water tank
- Charging station
- USB charging cable (USB-A to USB-C)
- 2 attachments
- 2 screws
- 2 wall plugs (Ø 5 mm)
- Operating instructions

Note

- Check the package for completeness and signs of visible damage.
- If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Customer Service hotline (see **20.1. Service** section).

7. Appliance description

(See fold-out page for illustrations)

- ❶ Nozzle
- ❷ Attachment
- ❸ Attachment socket
- ❹ Air supply openings
- ❺ Handpiece
- ❻ Water tank
- ❼ Filling opening, water tank/water tank lid
- ❽ Release button (for base cover)
- ❾ Base cover
- ❿ Hangers
- ⓫ Attachment holder
- ⓬ Recess for USB charging cable
- ⓭ USB-C port
- ⓮ Charging station
- ⓯ Charging contacts (charging station)
- ⓰ Charge indicator 
- ⓱ Indicator light  **soft**
- ⓲ Indicator light  **normal**
- ⓳ Indicator light  **jet**
- ⓴ Power level selection button 
- ⓵ On/Off button 
- ⓶ Release button
- ⓷ USB charging cable (USB-A to USB-C)
- ⓸ USB-C plug

8. Setting up the appliance / wall mounting

❶ Note

- The charging station ❶ may only be operated free-standing if used in dry rooms. Place the charging station ❶ on a horizontal, level and stable surface.
- In damp rooms, especially in the vicinity of water (bath), the charging station ❶ must be firmly fixed to the wall to protect it from falling into the water.
- Before wall mounting, make sure that there are no cables or pipes in the area where you want to hang the charging station ❶.

- 1) Mark two drill holes on the wall 3 cm apart and drill the marked holes with a 5 mm bit. Ensure that there is a mains power socket in the immediate vicinity.
- 2) Push the supplied wall plugs into the drill holes and screw the screws into the wall plugs. Take care that the screws still protrude about 3 mm out of the holes.
- 3) Feed the power cable through one of the recesses ❷ on the side of the charging station ❶, depending on which side the power socket is located.
- 4) Push the charging station ❶ onto the screws so that the screws protrude into the hangers ❿ at the back of the charging station ❶. Then press the charging station ❶ down so that the screws slide into the rails of the hangers ❿ and the charging station ❶ sits firmly on the screws.

9. Charging the appliance

⚠ ATTENTION!

- ▶ Due to the high current consumption, the appliance must be charged using a USB power adapter (not supplied). Do not charge the appliance from the USB port of a PC or notebook.
- ▶ To charge the appliance, use only a USB power adapter of protection class II approved for use with household devices (output voltage 5 V, max. 2 A).
- ▶ Use only the supplied USB charging cable ⑫ to charge the appliance.

ℹ Note

- ▶ The battery must be fully charged before using the appliance.
 - ▶ Before starting to use the appliance, check that it is in fault-free condition. You must not use the appliance if it has any faults.
 - ▶ Charge the appliance only when there is no water in the water tank ⑥.
 - ▶ A fully charged battery is enough for about 20 minutes of use.
- 1) Connect the USB-A plug on the USB charging cable ⑫ to a USB power adapter. Plug the USB power adapter into a power outlet.
 - 2) Connect the USB-C plug ⑭ of the USB charging cable ⑫ to the USB-C port ⑬ on the underside of the charging station ⑭.
 - 3) Place the handpiece ⑤ on the charging station ⑭ so that it slides onto the charging contacts ⑮ of the charging station ⑭.
 - 4) The charge indicator  ⑮ on the handpiece ⑤ lights up red during charging. When the battery is fully charged, the charge indicator  ⑮ lights up green.
 - 5) Place the handpiece ⑤ on the charging station ⑭ for charging after each use so that full power is always available.

10. Filling the water tank

⚠ ATTENTION!

- ▶ Fill the water tank ⑥ with pure tap water. If necessary, add a few drops of mouthwash. The water tank ⑥ should never be completely filled with additives such as mouthwash or mouth rinse.
 - ▶ Use only cold or lukewarm water with a maximum temperature of 40°C. Do not fill the water tank ⑥ with hot water.
- 1) Remove the handpiece ⑤ from the charging station ⑭. Hold the switched-off handpiece ⑤ horizontally.
 - 2) Open the water tank lid ⑦ and fill the tank with water through the filling opening ⑦.
 - 3) Close the water tank lid ⑦.

11. Fitting / removing attachments

ℹ Note

- ▶ Make certain that the handpiece ⑤ is switched off before fitting or removing the attachments ②.
- ▶ The appliance comes with two attachments ② with different coloured marking rings for different users.
- ▶ If used daily, replace the attachments ② after no more than six months.
- ▶ Use only the attachments ② recommended by the manufacturer. Information on ordering replacement attachments ② can be found in the section **21. Ordering replacement parts**. The attachment nozzles of Panasonic®* water jet flossers (nozzle model EW0950) can also be used for this appliance (*Panasonic is a registered trademark of Panasonic Corporation).

- ◆ To fit an attachment ② to the handpiece ⑤, insert it into the attachment socket ③. Press the attachment ② down until it completely engages in the attachment socket ③.
- ◆ To remove an attachment ② from the handpiece ⑤, press the release button ②② and pull the attachment ② upwards off the handpiece ⑤.
- ◆ Attachments ② that you are not currently using can be placed in one of the attachment holders ⑪ on the charging station ⑭.

12. Before first use

- 1) Rinse the attachments ② with clean water. You can place one attachment ② in the attachment holder ⑪ on the charging station ⑭. Fit the other attachment ② onto the handpiece ⑤.
- 2) Rinse the water tank ⑥ with clean water.
- 3) Fill the water tank ⑥ as described in the section **10. Filling the water tank**.
- 4) Point the nozzle ① of the attachment ② downwards into a washbasin. When doing so, make sure that you do not tilt the handpiece ⑤ too far forwards. Otherwise, the intake tube will no longer reach water or water could leak out of the air supply openings ④.
- 5) Press the On/Off button ① ②①. The handpiece ⑤ starts pumping water and discharging it via the nozzle ①. Let the handpiece ⑤ run until the water tank ⑥ is empty. This flushes all internal lines.
- 6) Press the On/Off button ① ②① to switch off the handpiece ⑤.

13. Power levels

This water jet flosser has three different power levels:

-  **jet**: powerful water jet, indicator light **jet** ①⑨ lights up.
For cleaning left over food from your interdental spaces.
-  **normal**: water jet mixed with air, indicator light **normal** ①⑩ lights up.
For cleaning left over food from your interdental spaces, for massaging your gums, for rinsing between teeth and gums.
-  **soft**: water jet mixed with a lot of air, indicator light **soft** ①⑪ lights up.
For gently massaging your gums.

Note

- If the appliance was completely discharged during the last use, it always starts on the power level **jet**.
- ◆ Press the power level selection button ① ②⑩ repeatedly until the desired power level is set.
- ◆ Press the On/Off button ① ②①. The appliance starts with the set power level.
- ◆ To switch off the water jet flosser, press the On/Off button ① ②① again.

14. Operation

CAUTION!

- The unfamiliar stimulation may cause some gum bleeding at first, which is not unusual. In this case, use a less intensive power level, such as **soft** or **normal**. If the bleeding gums last longer than two weeks, you should contact your dentist.
- Consult your dentist before using the appliance if you suffer from severe periodontitis, have injuries or abscesses in the mucous membrane of the oral cavity or underwent dental surgery within the last two months.

Note

- Use of the appliance is not a replacement for daily tooth cleaning with a toothbrush and toothpaste. Use the appliance only as a supplement to daily tooth cleaning with toothbrush and toothpaste.
- Make sure that you do not cover up the air supply openings **4** on the back of the handpiece **5** with your fingers during use.
- If you hold the handpiece **5** at an angle during use, water may leak out of the air supply openings **4**.
- The handpiece **5** switches off automatically after approx. 2 minutes at each power level, even if there is still water in the water tank **6**.

- 1) Fill the water tank **6**.
- 2) Place an attachment **2** on the handpiece **5**.
- 3) Press the power level selection button **20** to set the desired power level.
- 4) Place the attachment **2** in your mouth to prevent spray. Bend over the washbasin and hold the handpiece **5** as upright as possible (see Fig. 1). Direct the nozzle **1** of the attachment **2** towards your teeth/interdental spaces.

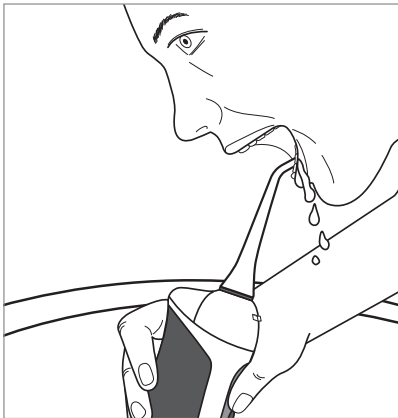
Fig. 1

- 5) Press the On/Off button **21** to start the handpiece **5**. Close your mouth a little to prevent water from spraying around.
- 6) To clean the interdental spaces, slowly guide the nozzle **1** of the attachment **2** from interdental space to interdental space on the inner and outer sides of your teeth. To massage your gums, point the nozzle **1** at the gums. Move the water jet slowly over your gums.
The length of the process will depend on how much cleaning the interdental spaces require. Continue cleaning until your teeth feel clean.
When massaging your gums, use the water jet flosser for as long as you find the massage comfortable.

Note

- The water in the water tank **6** is sufficient for approx. 45 seconds of use at power level **normal**.
 - If the water in the water tank **6** is used up in the meantime, stop the handpiece **5**. Remove it from your mouth and refill the water tank **6** through the filling opening **7**.
- 7) When the treatment is finished, press the On/Off button **21** to switch off the handpiece **5**.
 - 8) Only then should you take the attachment **2** out of your mouth and rinse it under clean, running water.

15. After use

Note

- Always empty the water tank **6** and all tubes after every use to prevent the build-up of bacteria.

When the handpiece **5** is switched off and you have finished the treatment:

- 1) Open the water tank lid **7** with the attachment **2** still fitted and pour out the rest of the water.
You can also remove the base cover **9** from the handpiece **5** to empty the water tank **6**. To do this, press the release button **8** and pull the base cover **9** off the handpiece **5**.
- 2) Close the water tank lid **7**.
- 3) Point the nozzle **1** of the attachment **2** downwards into a washbasin. Switch the handpiece **5** on for a couple of seconds to flush any remaining water from the water tank **6**.
- 4) Take the attachment **2** from the handpiece **5** and place it in the attachment holder **11** on the charging station **14**.
- 5) Place the handpiece **5** on the charging station **14**.

16. Cleaning and maintenance

WARNING!

-  Never immerse the handpiece **5** or charging station **14** in water or other liquids!

ATTENTION!

- Never open the housing of the charging station **14** or the handpiece **5** for cleaning!
- Do not use abrasive, aggressive or chemical cleaning agents. They could damage the outer surfaces of the appliance.
- ◆ Clean the charging station **14** only with a soft, lightly moistened cloth. Do not use detergents or solvents. They can damage the plastic surfaces.

- ◆ Remove the attachment **2** from the handpiece **5**. After use, clean the attachments **2** thoroughly under running water.
- ◆ Wipe off the handpiece **5** with a damp cloth. For more stubborn soiling, add a little detergent to the cloth and rinse it off afterwards with a cloth dampened with clean water. Dry the handpiece **5** afterwards.
- ◆ Remove the base cover **9** from the handpiece **5** by pressing the release button **8** and pulling the base cover **9** off the handpiece **5**. Rinse the water tank **6** thoroughly with warm water. Empty the water tank **6** completely into a sink after cleaning. If you do not use the appliance for more than a week, allow the inside of the water tank **6** to dry by laying the handpiece **5** flat without the base cover **9**. Allow the water tank **6** to air dry completely before refitting the base cover **9**.
- ◆ Replace the base cover **9** on the handpiece **5** and place it on the charging station **14**.

17. Storage

ATTENTION!

- Never store the device when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the device is stored for a longer period of time, the charge level of the battery must be checked regularly. The optimum charge level is between 50% and 80%.
- ◆ Disconnect the appliance from the power supply by disconnecting the USB charging cable **23** from the USB power adapter.
- ◆ Remove any water from the appliance and allow the appliance to dry completely (see section **16. Cleaning and maintenance**).
- ◆ Store the appliance in a dry and dust-free location out of direct sunlight.

18. Replacing the battery

CAUTION!

- ▶ Make sure that the battery is completely discharged before removing it. To do this, leave the appliance running until the battery is empty.

ATTENTION!

- ▶ This device contains a rechargeable battery that can only be replaced by trained personnel.

Note

- ▶ You need the following tools to change the battery:
 - 1 Phillips screwdriver size PH0
 - 1 flat-head screwdriver with a maximum blade thickness of 0.6 mm
 - 1 small side cutter
 - 1 soldering iron
- ▶ Ensure that the insulating rings are always fitted on the screws. Never insert the screws without the insulating rings. If you do, water may enter the appliance and damage it.
- ▶ Only use the batteries recommended by the manufacturer. Information on ordering replacement batteries can be found in the section **21. Ordering replacement parts**.

- 1) Remove the base cover  from the handpiece  by pressing the release button  and pulling the base cover  off the handpiece . Use the Phillips screwdriver to remove the two screws on the inside of the lower battery compartment cover. Then remove the lower battery compartment cover (see Fig. A).
- 2) Use the side cutter to cut the solder joints of the red and black wires on the positive and negative terminals of the lower battery compartment cover (see Fig. B). Make sure that you only cut the solder joints and not the wires.

- 3) Remove the top of the handpiece  by inserting the slotted screwdriver in the centre of the back of the handpiece  into the gap between the top and the handpiece . Release the top of the handpiece  with a little pressure and pull the top off the handpiece  (see Fig. C).
- 4) Use the Phillips screwdriver to remove the two lower screws on the upper battery compartment cover (see Fig. D). Slide the battery holder upwards out of the handpiece  (see Fig. E).
- 5) Disconnect the plug connection of the battery from the circuit board and then pull the battery out of the holder (see Fig. F).
- 6) Insert a new battery into the holder and connect the plug connection of the battery to the circuit board (see Fig. F).
- 7) Slide the battery holder from above into the handpiece  and use the Phillips screwdriver to tighten the two lower screws on the upper battery compartment cover.
- 8) Place the top on the handpiece  and press it firmly into the opening until you hear it click into place in the handpiece .
- 9) Use the soldering iron to solder the red wire to the positive terminal and the black wire to the negative terminal of the lower battery compartment cover. Pay attention to the labelling **+** and **-** on the lower battery compartment cover (see Fig. G).
- 10) Place the lower battery compartment cover on the handpiece  and use the Phillips screwdriver to tighten the two screws on the lower battery compartment cover.
- 11) Replace the base cover  on the handpiece  and place it on the charging station .

19. Disposal

19.1. Disposal of the appliance and packaging

Note

- ▶ Remove the internal battery before disposing of the appliance. To remove the battery, proceed as described in the section **18. Replacing the battery**. The battery must be disposed of in accordance with regulations. Ensure that the battery is disposed of in an officially designated place. Do not dispose of the battery in household waste.

The symbol of the crossed-out wheeled bin means that this appliance may not be disposed of in regular household waste at the end of its service life. The appliance must be deposited at an established collection point, recycling centre or disposal company.

Please erase all personal data before returning an appliance.

The packaging is made from environmentally friendly material which can be disposed of at your local recycling plant.

Dispose of the packaging in an environmentally friendly manner.

19.2. Disposal of batteries

Batteries/rechargeable batteries must be treated as hazardous waste and must therefore be disposed of in an environmentally sound manner by appropriate bodies (dealers, specialist dealers, public municipal bodies, commercial disposal companies). Batteries/rechargeable batteries may contain toxic heavy metals.

The heavy metals contained are identified by letters below the symbol: Cd = cadmium, Hg = mercury, Pb = lead.

For this reason, do not dispose of batteries/rechargeable batteries in domestic waste. Take them to a specialist collection point.

Only return batteries that are fully discharged.

20. Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 500431_2504 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 500431_2504.

20.1. Service

GB Service Great Britain
Tel.: 0800 404 7657
E-Mail: kompennass@lidl.co.uk

IE Service Ireland
Tel.: 1800 101010
E-Mail: kompennass@lidl.ie

MT Service Malta
Tel.: 80062230
E-Mail: kompennass@lidl.com.mt

IAN 500431_2504

20.2. Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

For EU market

KOMPERNASS HANDELS GMBH
BURGSTRASSE 21
44867 BOCHUM
GERMANY
www.kompennass.com

For GB market

Lidl Great Britain Ltd
Lidl House 14 Kingston Road
Surrey
KT5 9NU

21. Ordering replacement parts

You can order replacement parts for this product on the Internet at **www.kompennass.com**.

Scan the QR code with your smartphone/tablet.

You can use this QR

code to go directly to our website to view and order the available spare parts.

i NOTE

- If you have problems with your online order, you can contact our service centre by phone or e-mail.
- Always quote the article number (IAN) 500431_2504 in your order.
- Please note that online ordering of replacement parts is not possible for all countries.

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- Utilisez uniquement les embouts 2 recommandés par le fabricant. Vous trouverez des informations pour la commande d'embouts 2 de remplacement au chapitre **25. Commander des pièces de rechange**. Vous pouvez aussi utiliser pour cet appareil les embouts de remplacement des hydropulseurs Panasonic ©* (modèle EW0950) (*Panasonic est une marque déposée de Panasonic Corporation).