
www.lidl-service.com

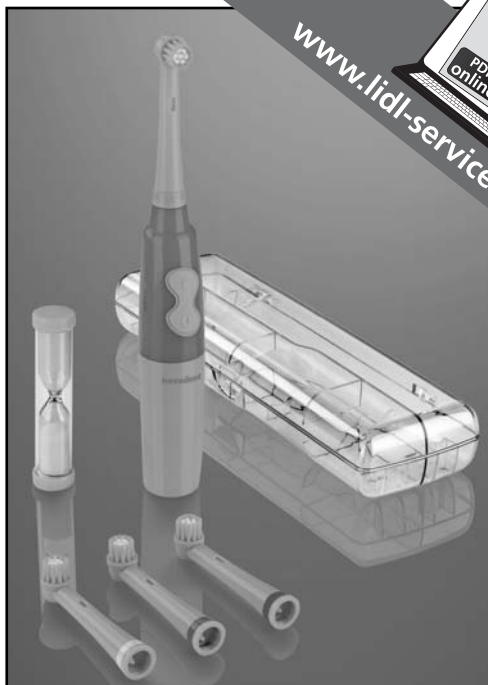KIDS' ELECTRIC TOOTHBRUSH NKZ 3 A1

(GB) (IE)

KIDS' ELECTRIC TOOTHBRUSH

Operating instructions

(FR) (BE)

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
(GB) (IE) (DK)
(NL) (BE)

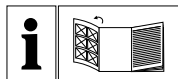GB IE

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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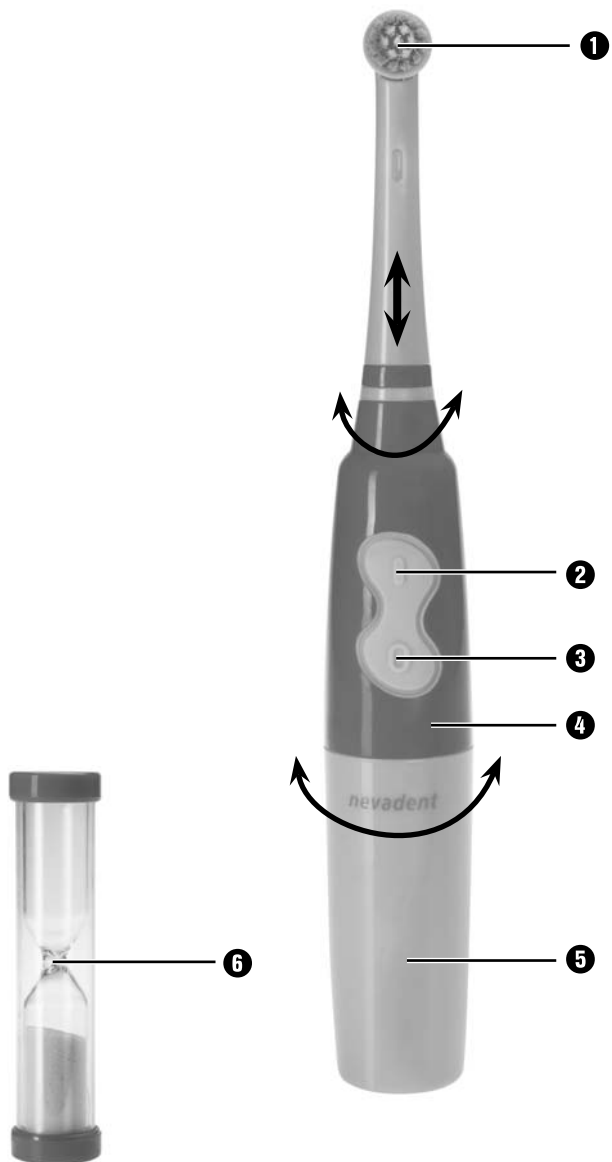
Contents

Introduction	2
Intended use	2
Package contents	2
Description of the appliance	2
Technical data	2
Safety instructions	3
Preparation	4
Fitting/removing brush heads	4
Brushing teeth	4
Cleaning and care	5
Disposal	6
Warranty and service	6
Importer	7

KIDS' ELECTRIC TOOTHBRUSH

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information on safety, usage and disposal. Before using the product, familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified areas of application. Please also pass on these operating instructions to any future owner.

Intended use

This appliance is intended exclusively for personal dental care at home. This appliance is intended solely for use in private households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Do not use it for commercial purposes.

Package contents

Hand element
4 brush heads
2 batteries
1 sticker sheet with 4 stickers
Sandglass
Storage box
Operating instructions

NOTE

Check contents for completeness and damage immediately after unpacking. Contact Customer Services if necessary.

Remove all packaging materials and plastic film from the appliance.

Description of the appliance

- ❶ Brush head
- ❷ On switch
- ❸ Off switch
- ❹ Hand element
- ❺ Battery compartment cap
- ❻ Sandglass

Technical data

Voltage supply	2 x 1.5 V 
Battery type	1.5 V, micro/AAA/LR03
Protection class	IPX4

Safety instructions

WARNING – RISK OF INJURY!

- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects penetrate the appliance.
- ▶ In the event of malfunctions or obvious damage, please contact our Customer Service department.
- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be performed by children unless they are under supervision.
- ▶ Children must not use the appliance as a plaything.

This appliance is not intended for children under 3 years of age.

- ▶ Do not allow the appliance to come into contact with hot surfaces. The appliance could be irreparably damaged!
- ▶ Do not open the appliance (except the battery compartment) yourself. Risk of injury if appliance is not properly assembled!

HANDLING BATTERIES

- ▶ Do not throw batteries into a fire. Do not recharge batteries. There is a risk of explosion and injury!
- ▶ Never open batteries and never solder or weld batteries. There is a risk of explosion and injury!
- ▶ Check the condition of the batteries at regular intervals. Leaking batteries can cause damage to the appliance.

HANDLING BATTERIES

- ▶ If you do not intend to use the appliance for an extended period, remove the batteries.
- ▶ If the batteries have leaked, use protective gloves. Clean the battery compartment and the battery contacts with a dry cloth.
- ▶ Keep batteries out of the reach of children. Children might put batteries into their mouth and swallow them. If a battery is swallowed, medical assistance must be sought immediately.

Preparation

Unscrew the battery compartment cap ⑤.

- 1) Insert the batteries as per the polarity diagram in the battery compartment.
- 2) Screw the battery compartment cap ⑤ back on.
- 3) Select a sticker and stick it on the front of the appliance.
- 4) Before initial use, clean the new brush heads ① with fresh water.

Fitting/removing brush heads

- Push a brush head ① onto the hand element ④ with a slight twisting motion.
- To remove the brush head ①, turn it a little anticlockwise while pulling it upwards and off.

NOTE

To avoid damage to the appliance, always use original brush heads.

(For replacement brush heads, see the order card)

Brushing teeth

⚠ WARNING!

Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.


- 1) Place the brush head ① onto the drive shaft on the handset ④.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

NOTE

To avoid spray, do not switch on the appliance until it is in your mouth

- 4) Turn the sandglass ⑥ over and press the on switch ②.

- 5) Divide your mouth into four parts (quadrants): upper right, upper left, lower right and lower left.

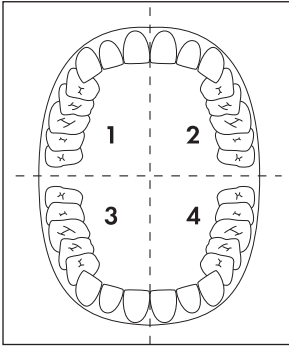- 6) Start in one quadrant:
Start with the brush head **1** on your gums, and move the brush upwards/downwards from the gums to the crown of the tooth. Proceed in this manner on both the inside and the outside of the tooth. Start again at the gums for each individual tooth.
Do not forget to brush the back of the rear molars.
To clean the chewing surfaces, guide the brush over them with a little pressure.
- 7) After about 30 seconds, move on to the next quadrant. Clean each jaw quadrant using the same technique.
- 8) After about two minutes, the sand in the sand-glass **6** will have run through, indicating that you have completed the two minutes of cleaning recommended by dental professionals.
- 9) Press the off switch **3**. The toothbrush stops.
- 10) Then take the toothbrush out of your mouth.

- 11) After use, rinse the brush head **1** thoroughly under running water.
- 12) Remove the brush head **1** from the handset **4**. To clean the drive shaft and the handset **4**, rinse them briefly with water. Hold the handset **4** so that the drive shaft is at a downward angle.

Cleaning and care

WARNING – PROPERTY DAMAGE!

-  Never submerge the hand element **4** in water or other liquids!
- Do not use detergents or solvents. They can damage the plastic surfaces.

- 1) Remove the brush head **1** from the handset **4**. After use, clean the brush heads **1** thoroughly under running water.
We recommend cleaning the brush heads **1** in a standard commercial denture cleaning fluid every 2–4 weeks. A new brush head **1** should be fitted after about 3 months, and the brush head must always be replaced when the bristles start to bend outwards.
- 2) To clean the drive shaft and the handset **4**, rinse them briefly with water. Hold the handset **4** so that the drive shaft is at a downwards angle.

Disposal

Never dispose of the appliance in the normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved waste disposal company or your municipal waste disposal facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.

Disposal of (rechargeable) batteries

Used batteries may not be disposed of in the household waste. Consumers are legally obliged to dispose of (rechargeable) batteries at a collection point in their community/city district or at a retail store. The purpose of this requirement is to ensure that batteries are disposed of in an environmentally friendly manner. Only dispose of batteries when they are fully discharged.

The packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

Warranty and service

This appliance is provided with a 3-year warranty valid from the date of purchase. The appliance has been manufactured with care and inspected meticulously prior to delivery. Please retain your receipt as proof of purchase. In the event of a warranty claim, please contact your Customer Service unit by telephone. This is the only way to guarantee free return of your goods.

The warranty only covers claims for material and manufacturing defects, not for damage sustained during carriage, components subject to wear and tear or damage to fragile components such as switches or batteries. This appliance is intended solely for private use and not for commercial purposes.

The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Your statutory rights are not restricted by this warranty.

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking, but no later than two days after the date of purchase. Repairs carried out after expiry of the warranty period shall be subject to charge.

GB Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: kompernass@lidl.co.uk

IAN 273033

IE Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: kompernass@lidl.ie

IAN 273033

Hotline availability:

Monday to Friday 08:00 - 20:00 (CET)

Importer

KOMPERNASS HANDELS GMBH

BURGSTRASSE 21

44867 BOCHUM

GERMANY

www.kompernass.com

ORDER CARD NKZ 3 A1

► www.kompernass.com

QUANTITY ORDERED

(max. 3 sets per order)

NAME OF ARTICLE

COST P. UNIT

TOTAL SUM

► Set of 8 Toothbrush heads

£ 5

£ _____
+ Postage and packaging £ 3

= £ _____

PAYMENT METHOD

☐ Bank Transfer

(Location, Date)

(Signature)

ORDER CARD NKZ 3 A1

► www.kompernass.com

1) Enter under "Sender/orderer" your name, address and telephone number (for possible queries) in block capitals.

2) Please transfer the complete amount to our account below:

Commerzbank AG, London

IBAN: GB96 COBA 4062 0130 4531 70

SWIFT (BIC): COBAGB2XXXX

Sort Code: 40-62-01

Account No: 30453170

Account Name: Kompernass Handelsgesellschaft mbH

■ On your transfer, please provide the article(s) you ordered as well as your complete name and address. Then post your completed order card to our postal address below.

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(Telephone)

IMPORTANT

- Please affix sufficient postage to the letter.
- Write your name as the sender on the envelope.

Order conveniently over the Internet
www.kompernass.com

ORDER CARD NKZ 3 A1

► www.kompernass.com

QUANTITY ORDERED
(max. 3 sets per order)

NAME OF ARTICLE

COST P. UNIT

TOTAL SUM

► Set of 8 Toothbrush heads

5,50 €

+ _____ €
Postage and packaging
3 €

= _____ €

PAYMENT METHOD

☐ Bank Credit Transfer

(Location, Date)

(Signature)

ORDER CARD NKZ 3 A1

✂ www.kompernass.com

1) Enter under "Sender/orderer" your name, address and telephone number (for possible queries) in block capitals.

2) Transfer the whole amount in advance to our account:

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SWIFT (BIC): PBNKDEFF 440

Bank: Postbank Dortmund AG

If you pay by bank credit transfer, on the remittance form please tick „shared“ (shared costs) in the box under the item `Costs`. Otherwise, we can not process your order.

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(Telephone)

IMPORTANT

► Please affix sufficient postage to the letter.

► Write your name as the sender on the envelope.

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