www.lidl-service.com

CEPILLO DENTAL ELÉCTRICO SPAZZOLINO DA DENTI ELETTRICO DAZD 3.7 Li B2

(ES)

CEPILLO DENTAL ELÉCTRICO

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(PT)

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(DE) (AT) (CH)

ELEKTRISCHE ZAHNBÜRSTE

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(IT) (MT)

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(GB) (MT)

ELECTRIC TOOTHBRUSH


Operating instructions

IAN 285254

(ES)

(IT)

(PT)

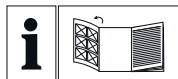ES

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Contents

Introduction	32
Proper use	32
Package contents	32
Appliance description	32
Technical details	32
Safety instructions	33
Set-up / wall bracket	34
First use	35
Fitting / removing brush heads	35
Operation	36
Brushing teeth	36
Travel lock	37
Cleaning and care	37
Travelling case	38
Disposal	38
Ordering spare parts	38
Kompernass Handels GmbH warranty	39
Service	40
Importer	40

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product.

The operating instructions are a constituent of this product. They contain important information about safety, use and disposal. Before using the product, familiarise yourself with all operating and safety instructions. Use this product only as described and only for the specified areas of application. Please also pass these operating instructions on to any future owner(s).

Proper use

This appliance is intended exclusively for personal dental care at home. This appliance is intended solely for use in domestic households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Not to be used for commercial purposes.

Package contents

Hand element
8 brush heads
Charger
2 screws
2 plugs
Travelling case
Operating instructions

NOTE

Check contents for completeness and damage immediately after unpacking. Contact Customer Services if necessary.

Remove all packaging materials and plastic film from the appliance.

Appliance description

- ❶ Brush head
- ❷ Lid
- ❸ On/Off switch
- ❹ "Clean" LED (orange)
- ❺ "Soft" LED (green)
- ❻ "Pulse" LED (yellow)
- ❼ Charge indicator /
Travel lock indicator 
- ❽ Charger
- ❾ Wall bracket

Technical details

Power supply:	100 - 240 V ~, 50/60 Hz
Power consumption:	2 W
Battery:	Li-ion rechargeable battery, 600 mAh
Protection rating (charger):	IP X7 Protection against temporary submersion
Protection class:	II / 
Charging time:	approx. 14 - 16 hours (initial charge)
Protection rating (handset):	IP X4 Protection against all-round splashes

Safety instructions

DANGER – RISK OF ELECTRICAL SHOCK!

- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects penetrate the appliance.
- ▶ In the event of malfunctions or obvious damage, please contact our Customer Service department.
- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way.
Risk of electric shock!
- ▶ The mains cable of this appliance cannot be replaced. In the event of the cable being damaged, the whole appliance must be scrapped.
- ▶ Never touch the charger if it has fallen into water. Remove the plug from the wall socket first! Risk of electric shock!

Do not submerge the handset or the charger in water or other liquids.

WARNING – RISK OF INJURY!

- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be performed by children unless they are under supervision.
- ▶ Children must not use the appliance as a plaything.

CAUTION – PROPERTY DAMAGE!


- ▶ The electric toothbrush may only be recharged with the supplied charger. The appliance could otherwise be irreparably damaged!
- ▶ Only operate the appliance with the original accessories supplied.

Set-up / wall bracket

NOTE

- ▶ The charger ❸ may only be operated free-standing if used in dry rooms.
- ▶ In damp rooms, especially in the vicinity of water (bath), the charger ❸ must be firmly fixed to the wall to protect it from falling into water.

- 1) Pull the wall bracket ❹ downwards from the charger ❸:

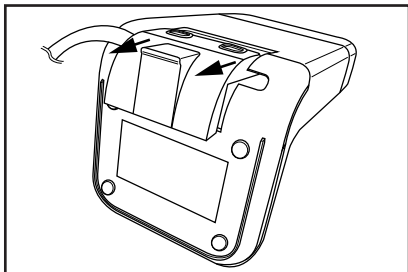- 2) Using the screw holes of the wall bracket ❹, mark the spacing for the boreholes on the wall. Check to ensure that there is a mains power socket in close vicinity.
- 3) Before drilling the holes, check to ensure that there are no electrical cables or water pipes behind the markings.
- 4) Use a 5 mm drill bit to drill the required holes and insert the wall plugs.
- 5) Affix the wall bracket ❹ using the screws provided.
- 6) After you have attached the wall bracket ❹, guide the power cable through the right or left cut-out in the charger ❸ (see magnification on the fold-out page). Using the guide rails, slide the charger ❸ onto the wall bracket ❹.
- 7) Ensure that the power cable is not trapped and that it lies properly in the cut-outs.

First use

NOTE

Before starting to use the appliance, check that it is in fault-free condition. If it is not, it must not be used.

- 1) Insert the plug into the mains power socket.
- 2) Before initial use, clean the new brush heads **1** with fresh water.
- 3) Open the lid **2** of the charger **8** and place the replacement brush heads on the brush holders. Re-close the lid **2**.
- 4) Place the handset on the charger **8** so that it slides onto the pin on the base.
- 5) Allow the handset to charge up for around 14 - 16 hours before using it for the first time. The charge indicator **7** on the front of the hand appliance lights up red during the charging process.
- 6) Use the handset repeatedly without charging the battery until the performance becomes noticeably weaker or the toothbrush stops working and the charge indicator **7** blinks red. Then recharge the battery for 14 - 16 hours. When the appliance is fully charged, the charge indicator **7** lights up green.
- 7) Charge and discharge the toothbrush completely 3 times. The battery will have then reached its full capacity.
- 8) You can now leave the handset in the charger **8** after every use, so as to always have full power available.

Fitting / removing brush heads

- When fitting the brush head **1**, attach it straight and with the bristles facing forwards towards the handset.
- To remove the brush head **1**, pull it upwards and off.

NOTE

For information on ordering replacement brush heads, read the section „Ordering replacement parts“.

You can also use the replacement brush heads from Oral-B®* electric toothbrushes for this appliance (except Pulsonic models).

*Oral-B is a registered trademark of Procter & Gamble Co.

Operation

This toothbrush has 3 different performing levels:

- “Clean”: normal speed,
orange “Clean” LED ④
For cleaning the teeth.
- “Soft”: slow speed,
green “Soft” LED ⑤
For cleaning the teeth of people with
sensitive gums.
- “Pulse”: alternating speeds,
yellow “Pulse” LED ⑥
For massaging the gums.

- Press the on/off switch ③ once to activate the
“Clean” mode. The “Clean” LED ④ lights up.
- Press once more to switch to “Soft” mode.
The “Soft” LED ⑤ lights up.
- Press it a third time to switch to “Pulse” mode.
The “Pulse” LED ⑥ lights up.
- To switch the toothbrush off, press the on/off
switch ③ once more.

Brushing teeth

⚠ WARNING!

Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.

- 1) Place the brush head ① on the handset.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

NOTE


To avoid spray, do not switch on the appliance until it is in your mouth

- 4) Use the On/Off switch ③ to select the desired operating mode (see section „Operation“).

NOTE

This toothbrush has a 30-second interval timer, so that all four jaw quadrants can be cleaned evenly.

- 5) Divide your mouth into four parts (quadrants):
upper left, upper right, lower left and lower right.

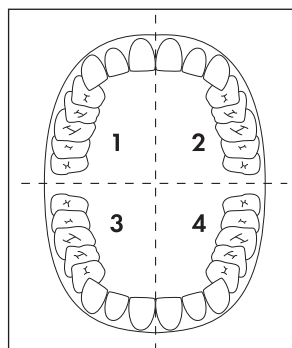- 6) Start in one quadrant:
Start with the brush head ① on your gums, and move the brush upwards/downwards from the gums to the crown of the tooth. Proceed in this manner on both the inside and the outside of the tooth. Start again at the gums for each individual tooth.
Do not forget to brush the back of the rear molars.
To clean the chewing surfaces, guide the oscillating brush head over the chewing surfaces with a little pressure.

- 7) After 30 seconds, the toothbrush stops briefly to signal that you should start cleaning the next quadrant. Clean each jaw quadrant using the same technique.
- 8) After two minutes the toothbrush stops its operation briefly 4 times. This informs you that you have complied with the dental recommendation of two minutes cleaning.
- 9) Press the On/Off switch ③ repeatedly until the toothbrush stops.
- 10) Then take the toothbrush out of your mouth.
- 11) After use, rinse the brush head ① thoroughly under running water.
- 12) Remove the brush head ① from the handset. To clean the drive shaft and the handset, allow a little water to flow over them. Hold the hand-set so that the drive shaft is pointing downwards.

NOTE

Regularly massage your gums. This promotes blood circulation and helps to keep them healthy. Set the mode to "Pulse" and guide the moistened brush over your gums with a light pressure.

Travel lock

When you store the toothbrush or wish to prevent it from switching on accidentally, you can activate the travel lock. The toothbrush cannot then be switched on by simply pressing the On/Off switch ③.

- To activate the travel lock, press and hold the On/Off switch ③ until the travel lock indicator ① ⑦ flashes green 5 times in a row. The On/Off switch ③ is now locked.

If you now attempt to switch on the toothbrush by briefly pressing the On/Off switch ③, the travel lock indicator ① ⑦ flashes green 5x to show that the travel lock is activated.

- To deactivate the travel lock, press and hold the On/Off switch ③ until the travel lock indicator ① ⑦ lights up 1x for a couple of seconds. The travel lock is now deactivated and you can start the appliance normally.

Cleaning and care

DANGER – ELECTRIC SHOCK!

- Always remove the plug of the charger ⑧ from the power socket before cleaning the appliance! Risk of electric shock!
-  Never submerge the handle or the charger ⑧ in water or other liquids.

NOTE

Never open the housing of the charger ⑧ or the toothbrush! They cannot be resealed in a water-tight manner.

- 1) Clean the charger ⑧ only with a soft, lightly moistened cloth. Do not use detergents or solvents. They can damage the plastic surfaces.
- 2) Remove the brush head ① from the handset. After use, clean the brush heads ① thoroughly under running water. We recommend cleaning the brush heads ① in a standard commercial denture cleaning fluid every 2 - 4 weeks. A new brush head ① should be fitted after about 3 months, and the brush head must always be replaced when the bristles start to bend outwards.

- 3) To clean the drive shaft and the handset, allow a little water to flow over them. Hold the handset so that the drive shaft is pointing downwards.

Travelling case

If you wish to take the toothbrush with you on a trip, you can pack both the hand element and two brush heads ❶ in the travelling case. When the battery is fully charged, you can use the toothbrush for 60 minutes without recharging. This is a guideline value only, depending on the pressure used when brushing your teeth, the maximum operating time may vary.

Disposal

Under no circumstance should the appliance be disposed of in the normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved disposal company or your municipal waste facility. Please observe the currently applicable regulations. Please contact your waste disposal centre if you are in any doubt.

NOTE

- The integrated rechargeable battery in this appliance cannot be removed for disposal.

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Ordering spare parts

You can order the following replacement parts for the product DAZD 3.7 Li B2:

- 8-pack of toothbrush heads

Order the replacement parts via the service hotline (see section "Service") or simply visit our website at www.kompernass.com:

NOTE

- Have the IAN number of the appliance ready (can be found on the cover of these operating instructions).

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on www.lidl-service.com.

Service

GB

Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: kompernass@lidl.co.uk

MT

Service Malta

Tel.: 80062230

E-Mail: kompernass@lidl.com.mt

IAN 285254

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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