SILVERCREST®

www.lidl-service.com

CELLULITE-MASSAGEGERÄT / CELLULITE MASSAGER APPAREIL ANTI-CELLULITE SMG 2.4 A1

(DE) (AT) (CH)

CELLULITE-MASSAGEGERÄT

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(FR) (BE)

APPAREIL ANTI-CELLULITE

Mode d'emploi

(PL)

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(SK)

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Contents

Introduction	20
Intended use	20
Package contents	21
Disposal of the packaging	22
Appliance description	23
Technical data	24
Safety information	25
Use	30
Cleaning	31
Storage	31
Disposal of the appliance	32
Mains adapter disposal	32
Kompernass Handels GmbH warranty	33
Service	36
Importer	36

Introduction

Congratulations on the purchase of your new appliance.
You have chosen a high-quality product.

The operating instructions are an integral part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

Intended use

This appliance is designed for massaging the skin on the buttocks, hips, upper thighs and upper arms. The massaging movements of the appliance promote the blood circulation in the corresponding skin areas.

This appliance is intended solely for personal, private use, it is not to be used for medical or commercial applications and is not a substitute for medical treatment. The appliance is designed for use in dry indoor rooms. It is not intended for use in the medical/therapeutic or commercial fields.

WARNING

Danger if not used as intended!

Failure to use the appliance for its intended purpose and/or any other type of use can be hazardous.

- ▶ Use the appliance only for its intended purpose.
- ▶ Observe the procedures described in these operating instructions.

Package contents

The appliance is supplied with the following components as standard:

- Cellulite massager
 - Mains adapter
 - Operating instructions
- 1) Remove all parts of the appliance and the operating instructions from the box.
 - 2) Remove all packaging material.

NOTE

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).

Disposal of the packaging

The packaging protects the appliance from damage during transport. The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable.

Returning the packaging into the material cycle conserves raw materials and reduces the amount of waste generated. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

- 1–7: Plastics,
- 20–22: Paper and cardboard,
- 80–98: Composites.

NOTE

- If possible, keep the original packaging during the warranty period so that the appliance can be properly packed for returning in the event of a warranty claim.

Appliance description

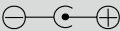
Figure A:

- ❶ Adjustable handle
- ❷ Switch: **0** = OFF, **1** = Slow rotation, **2** = Fast rotation
- ❸ Socket
- ❹ 6 removable massage rolls

Figure B:

- ❺ Plug
- ❻ Mains adapter

Technical data

Mains adapter ⑥	
Type	SW-0796
Input voltage	100 - 240 V ~ (alternating current), 50/60 Hz
Power consumption	0.7 A
Output voltage	12 V $\equiv$ (DC)
Output current	1.2 A
Protection class	II / $\square$ (double insulation)
Polarity	
Safety transformer	
Switched-mode power supply	
Protection type	IP20 Protection against solid objects with a diameter of more than 12.5 mm
Appliance	
Input voltage	12 V $\equiv$
Power consumption	1.2 A

Safety information

RISK OF ELECTRIC SHOCK

- ▶ Operate the appliance only with the supplied mains adapter.
- ▶ The appliance should only be connected to correctly installed and earthed mains power sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- ▶ Ensure that the power cable does not become wet or moist during use. Route the cable in such a way that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs to the appliance. All repairs must be performed by our Customer Service Unit or by a qualified specialist technician.

** RISK OF ELECTRIC SHOCK**

- ▶ Always unplug the appliance when it is not in use and before you clean it.
- ▶ Never lift, carry or drag the appliance by the power cable.
- ▶ To prevent accidents, have defective plugs, power cables or mains adapters replaced immediately by an authorised specialist technician or our Customer Service department.
- ▶ Use the appliance indoors only; never use it outside.
- ▶ Do not use the appliance while taking a bath or shower. Do not store the appliance anywhere where it could fall into the bath tub or the wash basin.
- ▶ Never try to retrieve an appliance that has fallen into water. Remove the mains adapter from the power socket first.

** RISK OF ELECTRIC SHOCK**

- ▶ Do not use the appliance if the power cable is damaged, the mains adapter is defective, the appliance is not working properly, if it has fallen into water or is otherwise damaged.
- ▶ Keep the power cable away from the rotating massage rollers.

** WARNING! RISK OF INJURY!**

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not allow children to use the appliance as a toy.
- ▶ People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- ▶ Do not use the appliance if you suffer from one or more of the following illnesses or complaints: circulatory problems, varicose veins, open wounds, bruising, skin tears, inflammation of the veins (phlebitis) or thrombosis.
- ▶ Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.
- ▶ Do not use the appliance under blankets or pillows.
- ▶ Do not use the appliance after taking substances that could reduce your capacity to react quickly (e.g. painkillers or alcohol).

⚠ WARNING! RISK OF INJURY!

- ▶ Stop treatment with the appliance immediately if you experience any discomfort or pain.
- ▶ Uninterrupted application of the appliance should not last any longer than 15 minutes. Afterwards, switch off the appliance for at least half an hour and let it cool down.
- ▶ Do not put your fingers or any other objects between the rotating massage rollers.
- ▶ Use the appliance only for its intended purpose.

CAUTION – PROPERTY DAMAGE!

Use the mains adapter/the appliance indoors only.

- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surface.

Use

WARNING – RISK OF INJURY!

- ▶ If you have long hair, tie it back before using the appliance on your upper arms. Otherwise, your hair could get caught in the appliance.
- ▶ Please beware of your long hair winded into the massage rolls **4** while using the machine.
- ◆ Insert the plug **5** on the mains adapter **6** into the socket **3** on the massager and then insert the mains plug **6** into a wall outlet.
- ◆ Slide your hand under the adjustable handle **1**. If this is too wide or too narrow, you can adjust it to fit your hand using the Velcro fastener.
- ◆ Press the massager gently against the part of the body that you want to massage.
- ◆ Switch on the appliance using the switch **2**. Select position **1** for a slow rotation of the rollers and position **2** for a fast rotation.

NOTE

- ▶ If you find the application unpleasant, try reducing the speed, adjust the pressure applied, change to a different part of the body or stop the massage completely.
- ◆ To switch the appliance off, push the switch **2** back to the **0** position.

Cleaning

RISK OF ELECTRIC SHOCK

- Before cleaning the appliance, always pull out the mains adapter **6** from the mains power socket.

CAUTION – PROPERTY DAMAGE!

- Do not use any aggressive cleaning agents or solvents. They could damage the surfaces.
- For easier cleaning, pull the massage rollers **4** off of the appliance. Clean them in warm water with a little mild soap. Rinse off the soap residues thoroughly and allow the massage rollers **4** to dry in the air. When they are dry, push them back onto the drive shaft of the massager.
- Wipe the appliance with a damp cloth. If necessary, add a little detergent to the cloth and wipe the appliance off afterwards with a cloth dampened with clear water.
- Wipe the mains adapter **6** and mains cable with a damp cloth. Make sure the mains adapter **6** has completely dried off before the next use.

Storage

- Clean the appliance as described in the section "Cleaning".
- Store the appliance in a clean and dust-free location.

Disposal of the appliance

Never dispose of the appliance in your normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance through an approved disposal facility or your municipal waste facility. Please observe the currently applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Mains adapter disposal

Do not dispose of the mains adapter in your normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

 Dispose of the mains adapter via an approved disposal centre or at your community waste facility. Please observe the currently applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Kompernass Handels GmbH

warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by telephone or by e-mail.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on www.lidl-service.com.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: kompernass@lidl.ie

IAN 311868

Importer

Please note that the following address is not the service address.
Please use the service address provided in the operating instructions.

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