# SILVERCREST®

[www.lidl-service.com](http://www.lidl-service.com)

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[www.kompernass.com](http://www.kompernass.com)

## Contents

|                                         |           |
|-----------------------------------------|-----------|
| <b>Introduction</b>                     | <b>34</b> |
| <b>Intended use</b>                     | <b>34</b> |
| <b>Package contents</b>                 | <b>35</b> |
| <b>Appliance description</b>            | <b>35</b> |
| <b>Safety information</b>               | <b>36</b> |
| <b>Function</b>                         | <b>39</b> |
| <b>Before first use</b>                 | <b>39</b> |
| <b>Use</b>                              | <b>40</b> |
| <b>Attachments</b>                      | <b>41</b> |
| <b>Programmes</b>                       | <b>41</b> |
| <b>Typical applications</b>             | <b>42</b> |
| Face                                    | 42        |
| Abdomen/waist                           | 42        |
| Abdomen/back                            | 43        |
| Legs                                    | 43        |
| Buttocks                                | 43        |
| Neckline                                | 44        |
| Arms                                    | 44        |
| <b>Cleaning</b>                         | <b>44</b> |
| <b>Storage</b>                          | <b>44</b> |
| <b>Disposal</b>                         | <b>45</b> |
| Disposal of the appliance/mains adapter | 45        |
| Disposal of the packaging               | 45        |
| <b>Appendix</b>                         | <b>46</b> |
| Technical data                          | 46        |
| Kompernass Handels GmbH warranty        | 47        |
| Service                                 | 48        |
| Importer                                | 48        |

## Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information on safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

## Intended use

This appliance is designed for vacuum massaging the skin.

The appliance is intended solely for personal, private use. It is not to be used for medical or commercial applications. Treatment with the appliance is no substitute for medical treatment. Do not use the appliance on animals.

The appliance is intended exclusively for use in dry indoor areas.

### **WARNING**

#### **Danger if not used as intended!**

Misuse of the appliance or use of the appliance for anything other than its intended purpose can be hazardous.

- ▶ Use the appliance only for its intended purpose.
- ▶ Observe the procedures described in these operating instructions.

## Package contents

The appliance is supplied with the following components as standard:

- Vacuum massager
- 2 attachments
- Mains adapter
- Storage bag
- Operating instructions

### DANGER

- ▶ Do not allow children to play with packaging materials.  
**There is a risk of suffocation.**

- ◆ Remove all parts of the appliance and the operating instructions from the box.
- ◆ Remove all packaging material.

### NOTE

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).

## Appliance description

(See fold-out page for illustrations)

- 1 Dial
- 2 On/Off switch
- 3 Programme switch
- 4 Control lamp
- 5 Connection for mains adapter
- 6 Storage bag
- 7 Mains adapter
- 8 Large attachment
- 9 Small attachment

## Safety information

### **RISK OF ELECTRIC SHOCK**

- ▶ Operate and charge the appliance only with the supplied mains adapter.
- ▶ The appliance should only be connected to correctly installed and earthed mains power sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- ▶ Ensure that the power cable does not become wet or moist during use. Route the cable in such a way that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not open the housing or attempt to repair the appliance yourself. This can lead to risks and also invalidates the warranty. Defective appliances should be repaired by authorised specialists only.
- ▶ In the event of a malfunction and before cleaning the appliance, disconnect the plug from the mains power socket.
- ▶ Never lift, carry or drag the appliance by the power cable.
- ▶ To prevent accidents, have defective plugs, power cables or mains adapters replaced immediately by an authorised specialist technician or our Customer Service department.
- ▶ Do not use the appliance while taking a bath or shower. Do not store the appliance anywhere where it could fall into the bath tub or the wash basin.

**⚠ RISK OF ELECTRIC SHOCK**

- ▶ Never try to retrieve an appliance that has fallen into water. Remove the mains adapter from the mains power socket first.
- ▶ Do not use the appliance if the power cable is damaged, the mains adapter is defective, the appliance is not working properly, if it has fallen into water or is otherwise damaged.
- ▶  Never immerse the appliance, the mains cable or the plug in water or any other liquids!

**⚠ WARNING! RISK OF INJURY!**

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- ▶ Do not allow children to use the appliance as a toy.
- ▶ Never use the appliance on animals.
- ▶ The appliance should not be used in the event of: serious cardiovascular disorders, chronic liver and kidney disease, skin diseases, varicose veins, open wounds, fresh operation scars, pregnancy.
- ▶ Do not use the appliance under blankets or pillows.
- ▶ Stop treatment with the appliance immediately if you experience any discomfort or pain.

## **⚠ WARNING! RISK OF INJURY!**

- ▶ Do not treat the same region of the body for too long or with too much force. Do not treat the same area several times in succession.
- ▶ Improper or excessive use may result in harmless but unattractive bruising.
- ▶ Use the appliance only for its intended purpose.

## **CAUTION – PROPERTY DAMAGE!**

- ▶  Only use the mains adapter and appliance indoors.
- ▶ Do not apply excessive amounts of cream or oil when using the appliance. This could adversely affect appliance functions and reduce the useful life.
- ▶ Do not use abrasive or caustic cleaning agents. These could damage the outer surfaces of the appliance.

## Function

Vacuum massage stimulates the tissue by the formation of negative pressure. The layers of skin are not pressed together as in conventional kneading massage, but instead pulled apart by the vacuum in the attachment.

In this way, deeper tissue layers can also be reached, while the blood and lymph circulation can be stimulated effectively. This helps to improve the nutrient supply to the skin as well as to detoxify the tissue.

## Before first use

- ◆ Clean the appliance as described in the section **Cleaning**.
- ◆ Charge the appliance fully before using it for the first time:
  - Connect the mains adapter **7** to the connection **5** of the appliance and plug the mains adapter **7** into a mains power socket. The control lamp **4** lights up red during charging.
  - Allow the appliance to charge for **four** hours and then disconnect it from the power supply if you do not wish to use it.

The appliance is now ready for use. Before starting, always read the **Use** section to learn how to handle the appliance safely and avoid injuries.

### NOTES

- ▶ The control lamp **4** goes out as soon as the appliance is fully charged. However, you should leave the appliance connected to the power supply for the recommended time before the first use.
- ▶ If the battery is discharged, the control lamp **4** will flash red. The appliance switches off after about 15 seconds. Plug the appliance into the power supply and recharge it.
- ▶ You can also use the appliance while it is charging. However, do not use it during the first charging process.

## Use

### **WARNING! RISK OF INJURY!**

To avoid bruises:

- ▶ Before switching on, make sure that appliance is set to the lowest intensity (dial **1** to "L" = LOW, low suction). Always start with the lowest intensity and increase this gradually if required (dial **1** direction "H" = HIGH, strong suction).
- ▶ Be especially careful when treating sensitive and non-firm regions of the body, such as the face!
- ▶ Do not treat the same region of the body for too long or with too much force. Do not treat the same area several times in succession.
- ▶ Every body reacts differently. If you find an application unpleasant, reduce the intensity, try another region of body or stop the treatment completely.

### **NOTE**

- ▶ Use the appliance every two days until you are satisfied with the result.
- ▶ Use the appliance only on regions of the body with sufficient fat tissue. Application on fat-free regions does not achieve any effect and will usually be experienced as unpleasant.
- ▶ Rub the areas to be treated with a small amount of a cosmetic product containing oil (e.g. body oil) so that the appliance slides over the skin more smoothly.
- ▶ You can stop the vacuum at any time by turning the dial **1**.

- ◆ Fit the appropriate attachment **8** **9**. See the **Attachments** section.
- ◆ Use the programme switch **3** to select the appropriate programme. See the **Programmes** section.
- ◆ Make sure the dial **1** is set to "L" (= LOW, low suction).
- ◆ Hold the opening of the attachment **8** **9** on the area to be treated and move the On/off switch **2** to "ON". The control lamp **4** lights up green in battery mode, or orange if the appliance is charging at the same time. The appliance immediately starts to suck in air and a vacuum is formed in the attachment **8** **9**.
- ◆ Move the appliance over your skin following the movements described in the **Typical applications** section.
- ◆ Change to an adjacent area in the suction pause. Make sure the areas do not overlap. You should not treat the same area several times in succession!

- ◆ To switch off the appliance, move the On/Off switch **2** to "OFF".  
The control lamp **4** goes out.
- ◆ Remove the attachment **8** **9** from the appliance.
- ◆ Clean the attachments used as described in the **Cleaning** section.

## Attachments

Two different attachments are included:

### Small attachment

The small attachment **9** is intended for treating small regions of the body, such as the face.

### Large attachment

The large attachment **8** is intended for treating larger areas such as the abdomen and back.

Proceed as follows to fit the attachments onto the appliance:

- ◆ Insert the attachment **8** **9** into the appliance head so that the locks on the attachment **8** **9** slide into the recesses in the appliance head and turn it clockwise until the attachment **8** **9** locks into place.
- ◆ To remove the attachment **8** **9** again, turn it anticlockwise until you can release it from the appliance head.

## Programmes

The appliance has four preset programmes, each with a different suction time:

### Programme 1 (P1)

2 seconds suction, 1 second suction pause

Select **P1** for all sensitive and non-firm regions of the body (e.g. face, insides of the upper arms etc.), especially for first-time treatments. If these areas respond well to the treatment, you can also use **P2** after a few attempts.

### Programme 2 (P2)

4 seconds suction, 2 seconds suction pause

Select **P2** if you are already familiar with the appliance and you think you can increase the suction time. If you are uncertain, it is better to remain with **P1**.

### Programme 3 (P3)

6 seconds suction, 2 seconds suction pause

Select **P3** for firm and tight regions of the body. If you are uncertain, it is better to remain with **P2**.

## Programme 4 (P4)

Continuous suction, no suction pause

This programme lets you adapt the treatment to your needs. However, only select **P4** if you are already experienced with the appliance and know how your body responds to the treatment. Turn the dial **1** to end the vacuum.

## Typical applications

This section illustrates various applications for possible parts of the body. Follow these examples for the best results.

Always start with a lower suction time (e.g. **P1** instead of **P2**) and the lowest suction intensity (dial **1** to "L"), especially if you are not very experienced with the appliance.

The illustrations show the treatment for just one half of the body. Repeat the treatment in exactly the same way for the other half of the body.

The arrows in the illustrations indicate the movement of direction you should follow for the treatment.

It may be advisable to ask a second person for help with areas that are difficult to reach, such as your back.

### Face

Select programme **P1** only, the facial area is usually very sensitive.


Always work from the centre of the face outwards.

### Abdomen/waist

Select programme **P1** or **P2** for the least firm regions.

Treat the remaining parts with **P2** or **P3**.

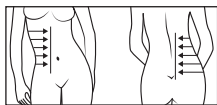Work horizontally from the outside to the centre of the body.

Work vertically from bottom to top.

## Abdomen/back

Select programme **P1** or **P2** for the least firm regions. Treat the remaining parts with **P2** or **P3**.

Start from the outside and work from bottom to top.

Start from the outside and work diagonally upwards towards the centre of the body. Finish off with an upward movement outwards.

## Legs

Select programme **P1** or **P2** for the least firm regions. Treat the remaining parts with **P2** or **P3**.

Work with a slightly diagonal movement downwards from outside to inside and finish off with an upwards movement.

Work horizontally from inside to outside and finish off with an upwards movement.

Start on the outside of the thigh and work diagonally upwards towards the centre of the body.


Start on the inside of the thigh and work in diagonal movements upwards towards the front and back of the leg.

## Buttocks

Select programme **P2** or **P3**.

Move the appliance in a curved motion from outside to inside on the relevant half of the buttocks.

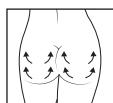Move the appliance in a curved motion from inside to outside on the relevant half of the buttocks.

## Neckline


Select programme **P1** or **P2** for the least firm regions.  
Treat the remaining parts with **P2** or **P3**.

Work from the bottom to the top.

## Arms

Select programme **P1** or **P2** for the least firm regions.  
Treat the remaining parts with **P2** or **P3**.

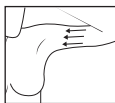Work from the elbow towards the armpit.

## Cleaning

### RISK OF ELECTRIC SHOCK

- ▶ Always disconnect the appliance from the mains adapter **7** before cleaning.

### CAUTION – PROPERTY DAMAGE!

- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surfaces.
- ◆ Always clean the attachments carefully with warm water and mild washing-up liquid after each use. Dry them thoroughly before storing or reusing them.
- ◆ Clean the appliance with a moist cloth. If required, add a little detergent to the cloth and wipe the appliance off afterwards with a cloth dampened with fresh water.
- ◆ Wipe the mains adapter **7** and mains cable with a damp cloth, if required. Make sure the mains adapter **7** and mains cable have completely dried off before the next use.

## Storage

- ◆ Clean the appliance as described in the section **Cleaning**.
- ◆ Put the appliance and accessories away in the storage bag provided **6** and store everything in a dry location.

# Disposal

## Disposal of the appliance/mains adapter

**Never dispose of the appliance or the mains adapter in your normal domestic waste.**

**These products are subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance and the mains adapter via an approved disposal company or your municipal waste facility. Please comply with all applicable regulations. Contact your waste disposal facility if you are in any doubt.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

### NOTE

- The integrated rechargeable battery in this appliance cannot be removed for disposal. To prevent risks, the removal or replacement of the battery may only be carried out by the manufacturer or the manufacturer's customer service department or a similarly qualified person. When disposing of the appliance, indicate that the appliance contains batteries.

## Disposal of the packaging

Returning the packaging to the material cycle conserves raw materials and reduces the amount of waste that is generated. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites

### NOTE

- If possible, keep the original packaging during the warranty period so that the appliance can be properly packed for returning in the event of a warranty claim.

Appendix

Technical data

|                                  |                                                                                                                                    |
|----------------------------------|------------------------------------------------------------------------------------------------------------------------------------|
| <b>Appliance</b>                 |                                                                                                                                    |
| Input voltage                    | 5 V  (direct current)                             |
| Power consumption                | 1 A                                                                                                                                |
| Polarity                         |                                                   |
| Rechargeable battery             | 3.7 V  (direct current)<br>800 mAh Li-ion battery |
| <b>Mains adapter</b>             |                                                                                                                                    |
| Type                             | DS006A0500100HE                                                                                                                    |
| Manufacturer                     | Dozen                                                                                                                              |
| Input voltage                    | 100–240 V ~ (alternating current),<br>50/60 Hz                                                                                     |
| Power consumption                | 0.15 A                                                                                                                             |
| Output voltage                   | 5.0 V  (direct current)                           |
| Output current                   | 1.0 A                                                                                                                              |
| Protection class                 | II /  (double insulation)                         |
| Efficiency class                 | 6                                                 |
| Nominal ambient temperature (ta) | 40 °C                                                                                                                              |
| Polarity                         |                                                 |
| Safety transformer               |                                                 |
| Switched-mode power supply       |                                                 |
| Protection type                  | IP 20: Protected against ingress of solid foreign objects with a diameter of > 12.5 mm.                                            |

## **Kompernass Handels GmbH warranty**

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

### **Warranty conditions**

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

### **Warranty period and statutory claims for defects**

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

### **Scope of the warranty**

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

## Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on [www.lidl-service.com](http://www.lidl-service.com).

## Service

### **GB** Service Great Britain

Tel.: 0800 404 7657

E-Mail: [kompernass@lidl.co.uk](mailto:kompernass@lidl.co.uk)

### **MT** Service Malta

Tel.: 80062230

E-Mail: [kompernass@lidl.com.mt](mailto:kompernass@lidl.com.mt)

IAN 313165

## Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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