MASSAGE GUN SMGB A1

GB IE NI

MASSAGE GUN

Short manual

DE AT CH

MASSAGE GUN

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IAN 391896_2201

GB IE NI



GB IE NI

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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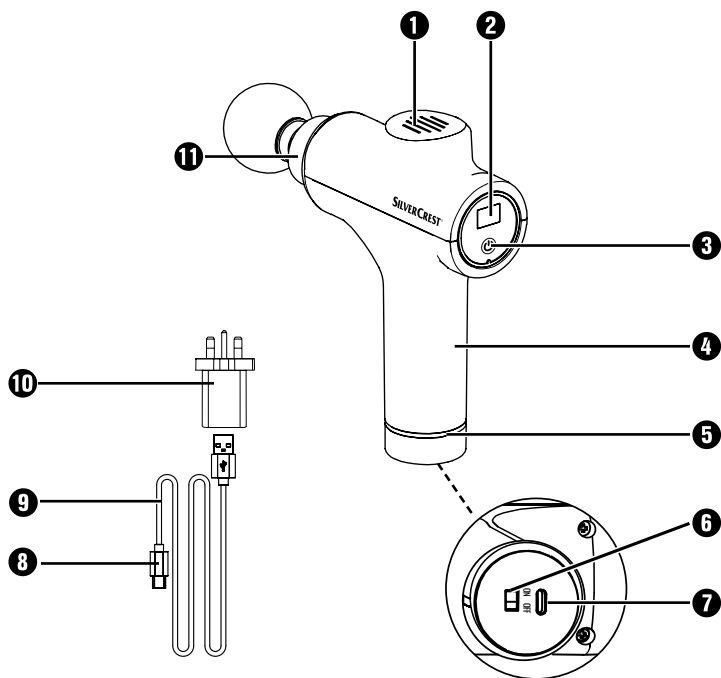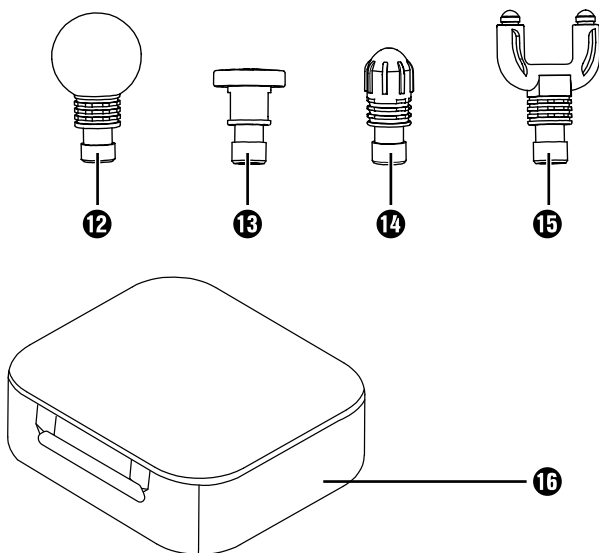A**B**

Contents

Information about this short manual	2
Intended use	2
Notes on trademarks	2
Safety information	3
Safety guidelines for rechargeable batteries	6
Package contents	7
Disposal of the packaging	7
Parts description	8
Before first use	8
Charging the appliance	8
Using the appliance	10
Selecting/fitting attachments	10
Examples of use	11
Switching the appliance on/off and selecting intensity	12
Intensity levels	13
Disposal	13
Technical specifications	14
Service	15
Importer	15

Information about this short manual

This document is a shortened print version of the complete operating instructions.

Scanning the QR code will take you directly to the Lidl Service page (www.lidl-service.com), where you can view and download the full version of the operating instructions by entering the article number (IAN) 391896_2201.

WARNING!

Observe the information in the complete operating instructions and the safety instructions to avoid personal injury and damage to property.

The short manual is an integral part of this product. Before using the product, please familiarise yourself with all operating and safety instructions. Keep the short manual in a safe place and hand over all documents when passing the product on to any future owner/user.

Intended use

This appliance and its attachments are intended exclusively for massaging the muscles of the neck, shoulders, arms, back (excluding the spine), thighs, calves and buttocks. Only use the attachments intended for this purpose. Never use the appliance directly on bones. The appliance is not intended for use on the face, head, neck or genital area.

This appliance is intended solely for personal, private use. It is not to be used for medical or commercial applications and is not a substitute for medical treatment. It is not intended for use in medical/therapeutic or commercial fields. Any other use or use beyond that specified is considered improper use.

Notes on trademarks

- USB® is a registered trademark of USB Implementers Forum, Inc.
All other names and products may be trademarks or registered trademarks of their respective owner

Safety information

To ensure safe operation of the appliance, follow the safety guidelines set out below.

RISK OF ELECTROCUTION!

- ▶ Charge the appliance only with the mains adapter supplied.
- ▶ Connect the mains adapter only to correctly installed and earthed mains sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- ▶ Make sure that the power cable does not become wet or moist. Route the cable in such a way that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs to the appliance. All repairs must be performed by our Customer Service Unit or by a qualified specialist technician.
- ▶ Never open the appliance housing. There are no parts inside the appliance which require user maintenance. In addition, this will invalidate the warranty.
- ▶ Switch off the appliance after each use and before cleaning it.
- ▶ Use the appliance in dry interior rooms only; never use it outside.
- ▶ Never immerse the appliance in water or other liquids!
- ▶ Do not use the appliance while taking a bath or shower. Do not store or keep the appliance in places from where it can fall into the bath tub, sink or other containers filled with water or other liquids.
- ▶ Never reach for the appliance if it has fallen into water.

RISK OF ELECTROCUTION!

- ▶ Do not use the appliance if the power cable is damaged, the mains adapter is defective, the appliance is not working properly, if the attachments show signs of damage, if it has fallen into water or is otherwise damaged.

Do not use the mains adapter if it has a damaged plug or cable. If the plug contacts of the mains adapter are defective or bent, the mains adapter cannot be repaired and must be scrapped.

WARNING! RISK OF INJURY!

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks must not be carried out by children unless they are supervised.
- ▶ Children must not be allowed to play with the appliance.
- ▶ Never leave children unattended with the packaging material. Keep children away from the packaging and the appliance. There is a risk of suffocation!
- ▶ People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- ▶ Do not use the appliance around your heart, face or genital area.
- ▶ Do not use the appliance if you suffer from one or more of the following diseases or medical conditions: circulatory disorders, varicose veins, open wounds, bruises, cracked skin, phlebitis, thrombosis, diabetes or slipped disc.

⚠ WARNING! RISK OF INJURY!

- ▶ Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.
- ▶ Do not use the appliance under blankets or pillows. Never use the appliance in the vicinity of petrol or other inflammable materials.
- ▶ Do not use the appliance while sleeping or after taking substances that may reduce your ability to react quickly (e.g. painkillers or alcohol).
- ▶ Never use the appliance directly on bones, e.g. on your spine, shins, the tops of your feet or similar places on your body.
- ▶ Use the appliance only on dry and clean parts of your body.
- ▶ Stop using the appliance immediately if you find the massage unpleasant or if you experience pain.
- ▶ Uninterrupted application of the appliance should not last any longer than 15 minutes. After use, switch off the appliance for at least 15 minutes and allow it to cool down.
- ▶ Do not insert your fingers or any other objects into the opening for the attachments. Keep long hair away from moving parts.
- ▶ Use the appliance only for its intended purpose.

CAUTION! PROPERTY DAMAGE!

Use the mains adapter/the appliance indoors only.

- ▶ Use the appliance only with the accessories supplied (mains adapter and attachments). Make sure that the attachments can always move freely.

CAUTION! PROPERTY DAMAGE!

- ▶ Protect the appliance and attachments from shocks, dust, chemicals or electromagnetic fields. Keep the appliance away from heat sources (radiators, ovens, etc.) and do not expose it to extreme fluctuations of temperature.
- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surface.
- ▶ Keep the ventilation slits free of dust and do not cover them.

NOTE

- ▶ No user action is required to switch the product between 50 and 60 Hz. The product switches automatically to either 50 or 60 Hz as required.

Safety guidelines for rechargeable batteries

The product is equipped or supplied with a lithium-ion battery. The battery in this appliance cannot be replaced.

⚠ RISK OF EXPLOSION!

- ▶ Do not open or short-circuit the battery. This could cause the battery to overheat and explode. There is a risk of fire!
- ▶ Never throw the battery into a fire or water.
- ▶ Do not subject the battery to mechanical stress.
- ▶ Charge the battery only with the supplied mains adapter (type GQ12-050200-ZB)
- ▶ Never leave the appliance unattended during charging.

⚠ WARNING! RISK OF INJURY DUE TO BATTERY LEAKAGE!

- ▶ High temperatures can affect the battery. Avoid direct sunlight and proximity to heat sources such as radiators.

⚠ WARNING! RISK OF INJURY DUE TO BATTERY LEAKAGE!

- ▶ Under extreme conditions, the battery cells may leak. Avoid contact with skin, eyes and mucous membranes. In the event of accidental contact, rinse the affected areas with plenty of clean water and contact a doctor.
- ▶ Leaking or damaged batteries can cause chemical burns to the skin. Always wear protective gloves if you have to touch the batteries.

Package contents

The appliance is supplied with the following components as standard:

- Massage gun
- 4 attachments (1 ball, 1 spiral, 1 U-shaped and 1 flat head)
- Storage bag
- Mains adapter with USB charging cable
- Short manual

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal, and are therefore recyclable.

Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

1–7: plastics, 20–22: paper and cardboard, 80–98: composites.

Parts description

Figure A (see the fold-out page):

- ❶ Ventilation slits
- ❷ Display/intensity display
- ❸ Button  (turn device on/off / set intensity level)
- ❹ Handle
- ❺ Ring display
- ❻ Main switch ON/OFF
- ❼ Connection socket
- ❽ USB plug (type C)
- ❾ USB charging cable (USB type-A to USB type-C)
- ❿ Mains adapter
- ⓫ Socket for attachments

Figure B:

- ❿ Ball head
- ⓫ Flat head
- ⓬ Spiral head
- ⓭ U-shaped head
- ⓮ Storage bag

Before first use

Before you can use the appliance, it needs to be charged for at least 3 ½ hours.

Charging the appliance

CAUTION

- Use only the supplied mains adapter ❿ to charge the battery.
- Do not use the appliance while it is charging.

- 1) Make sure that the appliance is switched off. The main switch ❻ must be set to OFF.

- 2) Insert the USB plug **8** into the USB connection socket **7** on the underside of the handle **4** (fig. 1):

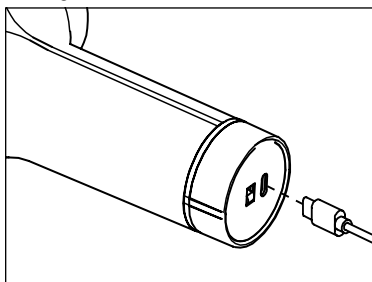
Fig. 1

- 3) Connect the USB charging cable **9** to the mains adapter **10** and plug the mains adapter **10** into a suitable mains socket.

✓ When the ring display **5** lights up green, the appliance is fully charged.

The battery level is displayed as follows:

Ring display 5	Meaning
Lights up red	Appliance switched on, battery charge less than 20 %.
Flashes red	Appliance is charging, battery level less than 30 %.
Flashes green	Appliance is charging, battery level at least 30 %.
Lights up green	Appliance switched on, battery fully charged.

- 4) Disconnect the appliance from the power supply and unplug the USB charging cable **9** from the appliance.

A full battery charge lasts up to 4 hours of operation at a medium intensity level. The battery can be recharged at any time and at any battery level.

CAUTION

- ▶ The appliance has a deep discharge protection device. When the battery charge level is too low, the display **2** shows the message **E 1**. Charge the appliance immediately to avoid damage to the battery.

Using the appliance

The massage gun provides a powerful, intensive massage for the thighs, calves, shoulders, arms and the back region (excluding the spine). Depending on the extent of the tension in the part of the body to be treated, a maximum treatment time of approx. 15 minutes is recommended.

A massage should always feel pleasant. A slight reddening of the treated areas and even slight bruising are a normal side effect, as the massage specifically promotes blood circulation. However, if you experience excessive skin reactions, discontinue use and consult a doctor.

Selecting/fitting attachments

- 1) Select the appropriate attachment for the desired treatment and body part as follows:

Attachment	Suitable for:
 12	Uniform and gentle massage to relax smaller muscle groups, e.g. upper arms, shoulders, waist and legs.
 13	Large areas all over the body, such as buttocks, thighs and back (excluding spine).
 14	Targeted and selective massage of deep muscle tissue, e.g. soles of the feet, palms of the hands, muscles of the shoulder blades.
 15	Uniform, 2-point stimulation of deep muscles in the lower back (except spine), neck, Achilles tendon.

- 2) Insert the desired attachment into the socket ⑪ and press it as far as it will go (Fig. 2):

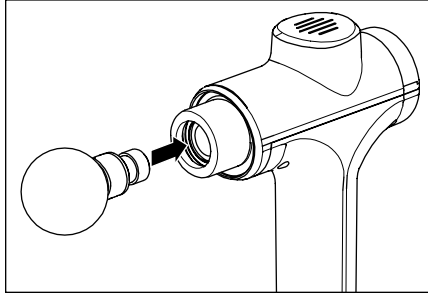
Fig. 2

- 3) To change the attachment, carefully pull it out of the socket ⑪ horizontally, turning it slightly if necessary.

Examples of use

NOTE

- Do not apply the attachments directly to a bone, only to muscle tissue.

Examples of application on arms (fig. 3), thighs (fig. 4), sole of the foot (fig. 5) and Achilles tendon (fig. 6).

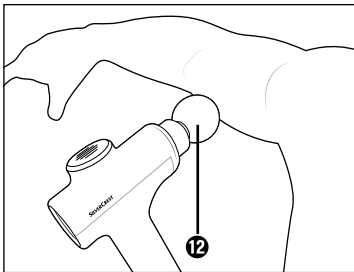
Fig. 3

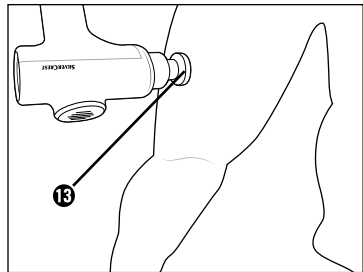
Fig. 4

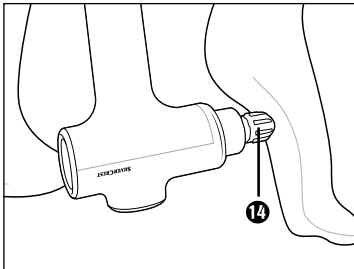
Fig. 5

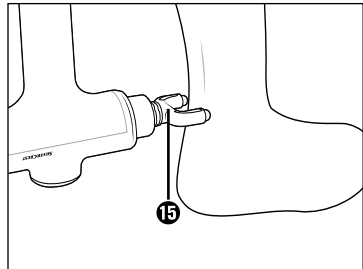
Fig. 6

Switching the appliance on/off and selecting intensity

- 1) Switch the appliance on by moving the main switch **6** on the underside of the handle **4** to the ON position.
 - ✓ The ring display **5** lights up green.
- 2) Press the button **1** **3** to start the appliance.
 - ✓ The display **2** shows the intensity level **1**.
The appliance starts on intensity level **1**.
- 3) To increase the intensity level, press the button **1** **3** briefly.
 - ✓ Each press increases the intensity by one level. This is shown on the display **2**.

NOTE

- After reaching the highest intensity level **05**, if you press the button **1** **3** again, the appliance stops working and switches into standby mode.
- 4) Massage the desired body part by gently moving the attachment back and forth over the muscle area.

NOTE

- Do not apply too much pressure, otherwise the appliance will automatically switch off to protect the motor from a possible overload. If this safety device has been triggered, **E3** will appear on the display **2**.
You need to restart the appliance.
- 5) After you have finished using the appliance, switch it off by pressing and holding the button **1** **3** for around 2 seconds and then move the main switch **6** to the OFF position.

NOTE


- The appliance is equipped with an automatic switch-off to prevent overheating. The appliance switches off automatically after 15 minutes of continuous operation. Allow the appliance to cool down for at least 15 minutes before using it again.

Intensity levels

The intensity of the massage varies depending on the level selected, as follows:

Level	Intensity
00	Off
01	1800 rpm - 30 Hz
02	2100 rpm - 35 Hz
03	2400 rpm - 40 Hz
04	2700 rpm - 45 Hz
05	3000 rpm - 50 Hz

Disposal

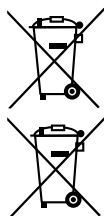The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If it is possible to do so without destroying the old appliance, remove the old batteries or rechargeable batteries before returning the appliance for disposal and take them to a separate collection point. In the case of permanently installed rechargeable batteries, you must indicate during disposal that the appliance contains a battery. The integrated rechargeable battery in this appliance cannot be removed for disposal.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Technical specifications

Mains adapter 10	
Manufacturers	DongGuan City GangQi Electronic Co.,Ltd. Privately Operated Industrial District, ShiYong, 532600 Heng Li Town, Dong Guan, Guangdong, P.R.CHINA Commercial registration number: 914419007929512928 Importer: KOMPERNASS HANDELS GMBH, BURGSTRASSE 21, 44867 BOCHUM, DEUTSCHLAND, Registergericht AG Bochum, Registernummer: HRB 4598
Model	GQ12-050200-ZB
Input voltage	100~240 V ~ (AC)
Input AC frequency	50/60 Hz
Output voltage	5.0 V $\equiv$ (DC)
Output current	2.0 A
Output power	10.0 W
Average operating efficiency	80.1%
No-load power consumption	0.06 W
Input current	max. 0.4 A
Efficiency class 6	
Polarity of the USB Port	P1  P4 
Safety transformer, short-circuit proof	
Switched-mode power supply	
Protection class	II /  (double insulation)
Nominal ambient temperature (ta)	45 °C
Protection type	IP20: Protection against ingress of solid objects with a diameter of more than 12.5 mm

Appliance

Input voltage	5.0 V $\equiv$ (DC)
Input current	2.0 A
Protection class	III /  (Protection by low voltage)
Integrated battery (Li-Ion)	2500 mAh / 7.4 V $\equiv$ (DC) / 18.5 Wh
	detachable supply unit

This mark confirms that the product complies with the product safety requirements applicable in the UK.

Service

Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

IAN 391896_2201

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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Lidl Great Britain Ltd