SHIATSU NECK MASSAGER SSMN 2 D2

(GB) (IE) (NI)

SHIATSU NECK MASSAGER

Short manual

(DE) (AT) (CH)

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IAN 391911_2201

(GB) (IE) (NI)



GB IE NI

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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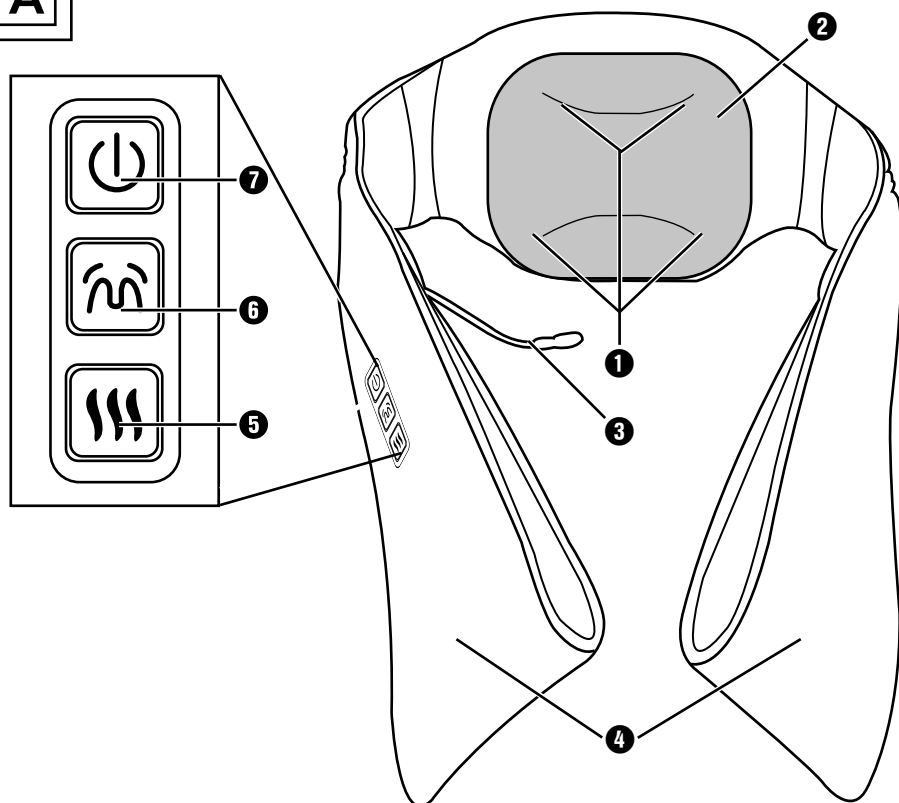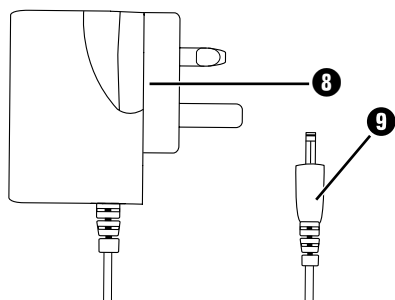A**B**

Index

Information about this short manual	2
Intended use	2
Items supplied	2
Disposal of the packaging	3
Appliance description	3
Technical specifications	4
Safety instructions	5
What is Shiatsu massage?	7
Operation	7
Connection	7
Fitting	7
Use	8
After the massage	9
Disposal of the appliance	9
Service	9
Importer	9

Information about this short manual

This document is a shortened print version of the complete operating instructions.

Scanning the QR code will take you directly to the Lidl Service page (www.lidl-service.com), where you can view and download the full version of the operating instructions by entering the article number (IAN) 391911_2201.

WARNING!

Observe the information in the complete operating instructions and the safety instructions to avoid personal injury and damage to property.

The short manual is an integral part of this product. Before using the product, please familiarise yourself with all operating and safety instructions. Keep the short manual in a safe place and hand over all documents when passing the product on to any future owner/user.

Intended use

This appliance is intended for massaging the neck and shoulders of people. This appliance is intended solely for personal use in the domestic sector, it is not to be used for medical or commercial applications and may not substitute for medical treatment. The appliance is designed for use in dry indoor rooms. It is not intended for use in the medical/therapeutic or commercial fields.

Items supplied

The appliance is delivered with the following components as standard:

- Shiatsu Neck Massager
- Mains adapter
- Short manual


NOTICE

- Check the contents to ensure that everything is present and for signs of visible damage.
- If the delivery is not complete or something is damaged due to defective packaging or through transportation contact the Service Hotline (see chapter **Service**).

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

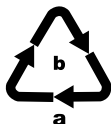Dispose of the packaging in an environmentally friendly manner.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

1–7: Plastics, 20–22: Paper and cardboard, 80–98: Composites.

Appliance description

See fold-out page for illustrations

Fig. A

- ➊ Massage heads
- ➋ Removable cover
- ➌ Socket
- ➍ Hand loops
- ➎ “Warmth” button  with LED light
- ➏ “Massage” button 
- ➐ “On/Off” button 

Fig. B

- ➑ Mains adapter
- ➒ Connection plug

Technical specifications

Mains adapter 	
Manufacturer:	Xiamen Keli Electronics Co., Ltd. No.17-19, XINGLIN NORTH 3RD ROAD, JIMEI DISTRICT, XIAMEN, FUJIAN PROVINCE, P.R. CHINA Commercial register number: 91350211737865674J Importer: KOMPERNASS HANDELS GMBH, BURGSTRASSE 21, 44867 BOCHUM, GERMANY, Register Court: AG Bochum, Register number: HRB 4598
Type	KL-AD-120100B
Input voltage	100-240 V ~ (AC)
Input AC frequency	50/60 Hz
Output voltage	12.0 V === (DC)
Output current	1.0 A
Output power	12.0 W
Average operating efficiency	84.1%
Low load efficiency (10%)	77.3%
No-load power consumption	0.08 W
Power consumption	0.7 A
Efficiency class 6	
Polarity	 (positive inside, negative outside)
Safety transformer, short-circuit proof	
Switched-mode power supply	
Protection class	II /  (double insulation)
Nominal ambient temperature (ta)	40 °C
Protection type	IP20 Protection against ingress of solid objects with a diameter of more than 12.5 mm

Appliance

Input voltage

12.0 V $\equiv$ (DC)

Power consumption

1.0 A

This mark confirms that the product complies with the product safety requirements applicable in the UK.

Safety instructions

RISK OF ELECTRIC SHOCK

- ▶ Operate the appliance only with the supplied power adapter KL-AD-120100B.
- ▶ Connect the appliance only to correctly installed and earthed mains power sockets. Ensure that the rating of the local power supply tallies completely with the details given on the rating plate of the appliance.
- ▶ Ensure that the power cable cannot become wet or damp while the appliance is in operation. Place it such that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs on the appliance. All repairs must be performed by our Customer Services or by qualified personnel.
- ▶ Always remove the plug from the mains power socket before you clean the appliance.
- ▶ Do not use the appliance in damp areas, such as in the bathroom.

Do not use the mains adapter if it has a damaged plug or cable. If the plug contacts of the mains adapter are defective or bent, the mains adapter cannot be repaired and must be scrapped.

- ▶ Do not insert needles or other sharp objects into the appliance.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not use an external timing switch or a separate remote control system to operate the appliance.
- ▶ This appliance may be used by children aged 8 over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use of the appliance safely and are aware of the potential risks. Children must not play with the appliance. Cleaning and user maintenance tasks must not be carried out by children unless they are supervised.
- ▶ Do not massage infants, young children, helpless persons or animals with the appliance.
- ▶ NEVER use the appliance if you have injuries or pain in the neck or cervical spine.
- ▶ Stop treatment with the appliance IMMEDIATELY if you feel it to be uncomfortable or painful.
- ▶ During use the appliance has a hot surface. Individuals who are insensitive to heat must take care when using the appliance.
- ▶ Do not use the appliance on swollen, burnt, inflamed, diseased or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the appliance.
- ▶ Ensure that your fingers or other body parts do not come between the rotating massage heads. Risk of crush injuries!
- ▶ The appliance may not be used in the face (e.g. eyes), on the throat, feet or other sensitive body parts.
- ▶ During use do not cover the appliance (with a pillow, blanket, etc.). NEVER use the appliance in the vicinity of petrol or other inflammable materials.

ATTENTION - PROPERTY DAMAGE!

Use the appliance only indoors.

NOTICE

- No user action is required to switch the product between 50 and 60 Hz. The product adapts automatically to either 50 or 60 Hz.

What is Shiatsu massage?

Shiatsu (shi = finger, atsu = pressure) is a form of massage developed in Japan, which itself emerged from traditional Chinese medicine. Shiatsu massage consists of a relaxing and powerful massage with fingers, thumbs or palms. This massager mimics the movements of shiatsu massage with the help of rotating massage heads ❶.

Operation

Connection

- 1) Connect the connection plug ❹ on the mains adapter ❸ with the socket ❸ on the appliance.
- 2) Connect the mains adapter ❸ into a mains power socket.

Fitting

- 1) Place the appliance around your neck so that the massage heads ❶ lie on the sides of your neck. You can use the hand loops ❷ to align it. The control panel with the buttons must be on your right.

Use

Ideally, you should lean back and relax during the massage. Provide for a relaxing atmosphere. For example, dim the lights a little.

NOTICE

- Ensure that you experience the massage as comfortable and relaxing. If needs be, vary the pressure and/or your position. Stop the massage immediately if it is painful or uncomfortable.
- You can also use the hand loops 4 to adjust the pressure, e.g. by pulling them a little tighter or leaving them looser.
- Never use the appliance for more than 15 minutes. A longer massage may lead to tension caused by overstimulation of the muscles.

You can select between 2 functions:

- Shiatsu massage with the massage heads ❶ (fast and slow)
- Warmth (2 of the massage heads ❶ illuminate and warm up)

You can combine various massage functions with each other or use them individually:

- Press the "On / Off" button  ❷. The massage heads ❶ begin with the slow massage.
- Press the button "Massage"  ❸ to increase the massage speed. Press the button "Massage"  ❸ once again to stop the massage heads ❶.
- Press the "Warmth" button  ❹ to switch the warmth on. The LED lights on the "Warmth" button  ❹ and the massage heads ❶ light up. Press the "Warmth" button  ❹ once again to switch the warmth function off. The lights on the massage heads ❶ and also the LED lights on the "Warmth" button  ❹ go out.

NOTICE

- If you have switched on the warmth function and the massage heads ❶ are rotating and you then stop the massage heads ❶ by pressing the "Massage" button  ❸, the warmth function also turns off at the same time. The LED lights on the "Warmth" button  ❹ flash. If you want to use the warmth function again, press the "Warmth" button  ❹ to reactivate the warmth function.

NOTICE

- The appliance switches off automatically after 15 minutes.
- If you want to terminate the massage before the end of 15 minutes and wish to switch the appliance off, press the "On/Off" button  ❷.

After the massage

- 1) Disconnect the connection plug ⑨ from the socket ③.
- 2) Remove the appliance.
- 3) Disconnect the mains adapter ⑧ from the mains power socket.

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Service

GB Service Great Britain
Tel.: 0800 404 7657
E-Mail: kompernass@lidl.co.uk

IE Service Ireland
Tel.: 1800 101010
E-Mail: kompernass@lidl.ie

IAN 391911_2201

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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Lidl Great Britain Ltd